


User Guide

Vaisala CARBOCAP[®] Carbon Dioxide Probe
GMP251 and GMP80P

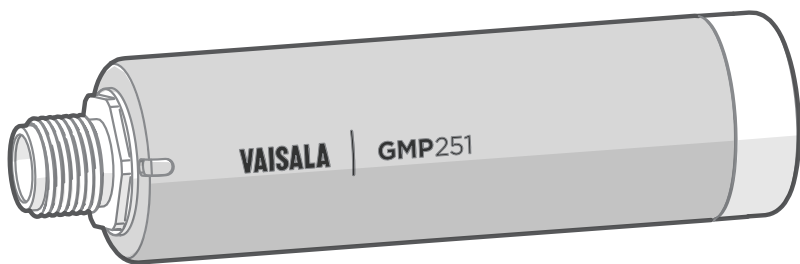

VAISALA

PUBLISHED BY

Vaisala Oyj
Vanha Nurmijärventie 21, FI-01670 Vantaa, Finland
P.O. Box 26, FI-00421 Helsinki, Finland
+358 9 8949 1
vaisala.com
docs.vaisala.com

© Vaisala 2024

No part of this document may be reproduced, published, or publicly displayed in any form or by any means, electronic or mechanical (including photocopying), nor may its contents be modified, translated, adapted, sold, or disclosed to a third party without prior written permission of the copyright holder. Translated documents and translated portions of multilingual documents are based on the original English versions. In ambiguous cases, the English versions are applicable, not the translations.

The contents of this document are subject to change without prior notice.

Local rules and regulations applicable to the products and services may vary, and they shall take precedence over the information contained in this document. Vaisala makes no representations on this document's compliance with the local rules and regulations applicable at any

given time, and hereby disclaims any and all responsibilities related thereto. You are instructed to confirm the applicability of the local rules and regulations and their effect on the intended use of the products and services.

This document does not create any legally binding obligations for Vaisala towards customers or end users. All legally binding obligations are set forth exclusively in the applicable contract or in the relevant set of General Conditions of Vaisala (vaisala.com/policies).

This product contains software developed by Vaisala or third parties. Use of the software is governed by license terms and conditions included in the applicable contract or, in the absence of separate license terms and conditions, by the General License Conditions of Vaisala Group.

Table of contents

1.	About this document.....	11
1.1	Version information.....	11
1.2	Related manuals.....	12
1.3	Documentation conventions.....	13
1.4	Trademarks.....	13
2.	Product overview.....	14
2.1	Introduction to GMP251.....	14
2.2	Basic features and options.....	15
2.2.1	Connectivity to Vaisala Insight software.....	15
2.2.2	Additional features with Indigo transmitters.....	16
2.2.3	Filter options.....	17
2.3	Operating principle of CO ₂ measurement.....	18
2.4	Environmental compensation.....	19
2.4.1	Temperature compensation.....	20
2.4.2	Pressure compensation.....	20
2.4.3	Background gas compensation.....	20
2.5	Probe startup.....	21
2.6	Filtering factor.....	21
2.7	Analog output overrange behavior.....	21
2.7.1	Analog output overrange example.....	22
2.8	Safety.....	23
2.8.1	ESD protection.....	24
3.	Installation.....	25
3.1	GMP251 probe dimensions.....	25
3.2	Recommended installation.....	25
3.3	Installation accessories.....	26
3.3.1	243261SP mounting flange.....	26
3.3.2	243257SP mounting clips.....	27
3.3.3	ASM213582 probe holder.....	28
3.3.4	ASM212017SP spray shield.....	29
3.4	Power supply.....	29
3.5	Wiring.....	30
4.	Modbus.....	31
5.	Vaisala Insight PC software.....	32
5.1	Connecting to Insight software.....	32
6.	Using GMP251 with Indigo200 series transmitters.....	34
6.1	Introduction to Indigo200 transmitters.....	34
6.2	Connecting probes.....	35
6.3	Analog output configuration overview.....	36
6.3.1	Analog output configuration example.....	37

6.4	Modbus overview.....	37
6.4.1	Default communication settings.....	39
6.5	Relay configuration overview.....	39
6.6	Configuring filtering factor with Insight.....	41
7.	Using GMP251 with Indigo80 handheld indicator.....	43
7.1	Indigo80 handheld indicator.....	43
7.1.1	Probe compatibility.....	44
7.1.2	Indigo80 keypad.....	44
7.1.3	Indigo80 main menu.....	44
7.2	Connecting probes to Indigo80.....	45
7.3	Configuring probe features with Indigo80 handheld indicator.....	47
7.3.1	Example of configuring filtering factor with Indigo80 handheld indicator.....	47
7.4	Configuring environment settings with Indigo80 handheld indicator....	48
8.	GMP80P Portable CO₂ Probe with Pump Sampling.....	50
8.1	GMP80P product overview.....	50
8.2	GMP80P probe structure.....	51
8.3	Connecting GMP80P to Indigo80.....	51
8.3.1	Configuring GMP80P for sampling use.....	52
8.4	GMP80P installation and sampling.....	53
8.4.1	Nafion membrane tube in sampling.....	53
8.4.2	CO ₂ dilution coefficients for pump sampling method.....	54
8.4.3	Sampling using plastic tube.....	57
8.4.4	Sampling using steel pipe.....	58
9.	Operating with MI70 indicator.....	59
9.1	Overview of MI70 support.....	59
9.1.1	Probe warm-up and MI70 automatic power off.....	59
9.2	Basic display.....	59
9.3	Graphical display.....	59
9.4	Main menu.....	60
9.5	Connecting probe to MI70 indicator.....	60
9.6	MI70 indicator parts.....	61
9.7	Holding and saving the display.....	61
9.8	Recording data.....	62
9.9	Changing environmental compensation settings with MI70 indicator....	62
10.	Vaisala Industrial Protocol.....	64
10.1	Vaisala Industrial Protocol overview.....	64
10.2	Serial interface settings.....	64
10.3	Physical interface.....	64
10.4	Connecting with a computer.....	64
10.4.1	Installing driver for the USB service cable.....	66
10.5	Accessing serial commands from Modbus or analog mode.....	67
10.6	Enabling Modbus mode from Vaisala Industrial Protocol.....	67
10.7	Changing from digital output to analog output.....	68

10.8	Serial commands.....	68
10.9	Device information and status.....	70
10.10	Serial line output and communication.....	73
10.11	Analog output.....	81
10.12	Calibration and adjustment commands.....	84
10.13	Environmental compensation commands.....	88
10.14	Other commands.....	93
11.	Calibration and adjustment.....	94
11.1	Calibration setup.....	94
11.2	Effect of environmental compensations.....	95
11.3	Limits of adjustment.....	95
11.4	Adjustment types.....	95
11.5	DRW244827SP calibration adapter.....	97
11.6	Overview of calibration with Indigo200 using Insight.....	98
11.6.1	Environmental compensation.....	99
11.6.2	CO ₂ probe adjustment example.....	100
11.7	Calibration and adjustment with Indigo80 handheld indicator.....	101
11.7.1	Performing guided calibration with Indigo80.....	102
11.7.2	Performing manual adjustment with Indigo80 using predetermined values.....	103
11.7.3	Restoring factory adjustment with Indigo80 handheld indicator..	104
11.8	Calibration and adjustment with MI70 indicator.....	104
11.8.1	1-point adjustment with an MI70-compatible reference probe.....	105
11.8.2	1-point adjustment with a reference gas.....	106
12.	Maintenance.....	109
12.1	Cleaning GMP251 and GMP80P.....	109
12.1.1	Chemical tolerance.....	109
12.2	Detaching GMP251 from GMP80P.....	109
12.3	Changing the GMP251 filter.....	110
13.	Troubleshooting.....	111
13.1	Problem situations with GMP251.....	111
13.2	Error messages.....	111
13.3	Analog output error state.....	113
14.	Technical data.....	114
14.1	Product specifications.....	114
14.1.1	GMP251 specifications.....	114
14.1.2	GMP80P specifications.....	117
14.2	Spare parts and accessories.....	120
14.2.1	GMP251 spare parts and accessories.....	120
14.2.2	GMP80P spare parts and accessories.....	121

14.3	Dimensions.....	121
14.3.1	GMP251 probe dimensions.....	121
14.3.2	GMP80P dimensions.....	122
14.3.3	243261SP mounting flange dimensions.....	122
14.3.4	DRW244827SP calibration adapter dimensions.....	123
14.3.5	ASM212017SP spray shield dimensions.....	123
Appendix A: Modbus reference.....		125
A.1	Default communication settings.....	125
A.2	Function codes.....	125
A.3	Data encoding.....	125
A.3.1	32-bit floating point or 32-bit integer format.....	125
A.3.2	16-bit integer format.....	126
A.4	Modbus registers.....	126
A.4.1	Measurement data registers.....	127
A.4.2	Configuration registers.....	127
A.4.3	Status registers.....	129
A.5	Device identification objects.....	130
A.6	Exception responses.....	130
A.7	Modbus communication examples.....	131
A.8	Filtering factor.....	133
Appendix B: CO₂ dilution coefficients reference.....		134
B.1	CO ₂ dilution coefficients.....	134
Warranty.....		135
Technical support.....		135
Recycling.....		135

List of figures

Figure 1	GMP251 probe parts.....	14
Figure 2	Probe cuvette with mirror and sensor chips.....	18
Figure 3	CO ₂ measurement in the measurement cuvette.....	19
Figure 4	Example of analog output overrange behavior.....	23
Figure 5	GMP251 dimensions.....	25
Figure 6	Probe with 243261SP mounting flange.....	26
Figure 7	Probe in 243257SP mounting clips.....	27
Figure 8	Probe holder ASM213582 dimensions.....	28
Figure 9	Probe with ASM212017SP spray shield.....	29
Figure 10	Connecting probe to Insight.....	32
Figure 11	Indigo201 with and without display.....	34
Figure 12	Connecting probes and cables to Indigo200.....	35
Figure 13	Analog output configuration options example in Insight.....	36
Figure 14	Analog output type (current or voltage) selection in Insight.....	37
Figure 15	Insight PC software Communication menu view.....	39
Figure 16	Relay configuration options (CO ₂ probe example).....	40
Figure 17	Relay icons on display (relay A active, relay B not active).....	41
Figure 18	Filtering factor setting in Insight.....	41
Figure 19	Indigo80 handheld indicator.....	43
Figure 20	Indigo80 keypad.....	44
Figure 21	Indigo80 main menu.....	44
Figure 22	Connecting probes to Indigo80.....	45
Figure 23	Vaisala GMP80P CO ₂ probe.....	50
Figure 24	GMP80P probe structure.....	51
Figure 25	GMP80P connected to Indigo80.....	52
Figure 26	Nafion tube, Vaisala item code 212807SP.....	54
Figure 27	CO ₂ sampling from incubator using tube, with GMP80P connected to Indigo80.....	57
Figure 28	CO ₂ sampling from incubator using pipe.....	58
Figure 29	MI70 basic display.....	59
Figure 30	MI70 indicator parts.....	61
Figure 31	CO ₂ reading with Tcomp and Pcomp on MI70 screen.....	62
Figure 32	Probe compensation settings on MI70 screen.....	63
Figure 33	Example of analog output overrange behavior.....	83
Figure 34	DRW244827SP calibration adapter with probe inserted.....	97
Figure 35	Insight calibration menu main view (CO ₂ probe).....	98
Figure 36	Insight Calibrate menu selection.....	99
Figure 37	Insight Restore factory adjustment selection.....	99
Figure 38	Insight environmental compensation configuration example with CO ₂ probe.....	100
Figure 39	Calibration methods available with Indigo80.....	101
Figure 40	Opening the filter.....	110
Figure 41	GMP251 dimensions.....	121
Figure 42	GMP80P dimensions, front and side.....	122
Figure 43	243261SP mounting flange dimensions.....	122

Figure 44 243261SP mounting flange dimensions, cross section..... 122

Figure 45 DRW244827SP calibration adapter dimensions..... 123

Figure 46 ASM212017SP spray shield dimensions..... 123

Figure 47 GMP251 spray shield cross section..... 124

List of tables

Table 1	Document versions (English).....	11
Table 2	Related manuals.....	12
Table 3	Analog output overrange clipping and error limits.....	22
Table 4	M12 male connector.....	30
Table 5	Default Modbus serial communication settings.....	31
Table 6	Indigo202 Modbus registers.....	38
Table 7	Default Modbus serial communication settings.....	39
Table 8	CO ₂ dilution coefficients for measurements with Nafion tube, T _d of sampled gas 0 ... +60 °C.....	55
Table 9	CO ₂ dilution coefficients for GMP80P sampling methods at 25 °C (77 °F) ambient temperature.....	56
Table 10	Default serial interface settings.....	64
Table 11	Basic serial commands.....	68
Table 12	Advanced serial commands.....	69
Table 13	? command.....	70
Table 14	errs command.....	71
Table 15	help command.....	71
Table 16	snum command.....	71
Table 17	system command.....	72
Table 18	time command.....	72
Table 19	vers command.....	72
Table 20	addr command.....	73
Table 21	close command.....	73
Table 22	form command.....	73
Table 23	Output parameters for form command.....	75
Table 24	Modifiers for form command.....	75
Table 25	intv command.....	76
Table 26	open command.....	76
Table 27	r command.....	76
Table 28	s command.....	77
Table 29	sdelay command.....	77
Table 30	send command.....	77
Table 31	seri command.....	78
Table 32	smode command.....	79
Table 33	amode command.....	81
Table 34	aover command.....	82
Table 35	asel command.....	83
Table 36	adate command.....	84
Table 37	atext command.....	85
Table 38	cdate command.....	85
Table 39	ctext command.....	85
Table 40	cco2 command.....	86
Table 41	env command.....	88
Table 42	o2cmode command.....	90
Table 43	pcmode command.....	91

Table 44	rhcmode command.....	91
Table 45	tcmode command.....	92
Table 46	frestore command.....	93
Table 47	pass command.....	93
Table 48	reset command.....	93
Table 49	GMP251 measurement performance.....	114
Table 50	GMP251 operating environment.....	115
Table 51	GMP251 inputs and outputs.....	116
Table 52	GMP251 mechanical specifications.....	116
Table 53	GMP251 compliance.....	117
Table 54	GMP80P measurement performance.....	117
Table 55	GMP80P operating environment.....	118
Table 56	GMP80P powering.....	118
Table 57	GMP80P digital communication.....	118
Table 58	GMP80P mechanical specifications.....	118
Table 59	GMP80P compliance.....	119
Table 60	GMP80P manufacturer contact information.....	119
Table 61	GMP251 spare parts and accessories.....	120
Table 62	GMP80P spare parts and accessories.....	121
Table 63	Default Modbus serial communication settings.....	125
Table 64	Modbus function codes.....	125
Table 65	Interpretation of 16-bit signed integer values.....	126
Table 66	Modbus measurement data registers (read-only).....	127
Table 67	Modbus configuration data registers (writable).....	127
Table 68	Modbus status registers (read-only).....	129
Table 69	Device identification objects.....	130
Table 70	Modbus exception responses.....	130
Table 71	CO ₂ dilution coefficients for measurements with Nafion tube, T _d of sampled gas -40 ... +60 °C.....	134

1. About this document

1.1 Version information

Table 1 Document versions (English)

Document code	Date	Description
M211799EN-J	May 2024	<p>Added sections:</p> <ul style="list-style-type: none"> • GMP80P Portable CO₂ Probe with Pump Sampling (page 50) • Performing guided calibration with Indigo80 (page 102) • Detaching GMP251 from GMP80P (page 109) • GMP80P specifications (page 117) • GMP80P spare parts and accessories (page 121) • GMP80P dimensions (page 122) <p>Updated sections:</p> <ul style="list-style-type: none"> • Trademarks (page 13) • Basic features and options (page 15) • Additional features with Indigo transmitters (page 16) • Configuring environment settings with Indigo80 handheld indicator (page 48) • Calibration setup (page 94) • Performing manual adjustment with Indigo80 using predetermined values (page 103) • Cleaning GMP251 and GMP80P (page 109) • Status registers (page 129): removed register 2052 and the related table <i>Error codes in register 0803hex</i>. • CO₂ dilution coefficients (page 134) • Updated Modbus register number notation throughout the document
M211799EN-H	July 2023	<ul style="list-style-type: none"> • Added information on using the GMP251 probe with the Indigo80 handheld indicator. • References to Indigo200 transmitters' wireless configuration interface replaced with information about the new Indigo200 transmitters that have a USB-C interface. For instructions on configuring GMP251 probes with the earlier wireless configuration interface, see <i>GMP251 User Guide</i> revision M211799EN-G, available at docs.vaisala.com. • Added section <i>ASM213582 probe holder</i> (page 28). • Updated sections: <ul style="list-style-type: none"> • Recommended installation (page 25) • Environmental compensation commands (page 88) • Error messages (page 111) • GMP251 specifications (page 114) • GMP251 spare parts and accessories (page 120) • Status registers (page 129): added register 2052. • Removed section <i>Patent notice</i>.

Document code	Date	Description
M211799EN-G	March 2018	<ul style="list-style-type: none"> Added information about using the probe with Vaisala Insight PC software. PuTTY terminal settings instructions updated. Calibration uncertainty specification updated. Added recommendation to disable MI70 indicator's automatic power off feature when used with GMP251. Updated instructions on using MI70 to configure environmental compensations. Default device address in serial line communication corrected (240). Corrected the unit information of sdelay command. Added a note about maximum voltage output to analog error state output configuration. Information on Modbus register numbers updated.

1.2 Related manuals



For the latest versions of these documents, see docs.vaisala.com.

Table 2 Related manuals

Document code	Description
M211798EN	Vaisala GMP251 Carbon Dioxide Probe Quick Guide
M211897EN	Vaisala GMP252 Carbon Dioxide Probe User Guide
M211893EN	Vaisala GMP252 Carbon Dioxide Probe Quick Guide
M212945EN	Vaisala GMP80P Portable Carbon Dioxide Probe with Pump Sampling Quick Guide
M212722EN	Vaisala Indigo80 Handheld Indicator User Guide
M212842EN	Vaisala Indigo201 Analog Output Transmitter User Guide
M212815EN	Vaisala Indigo201 Analog Output Transmitter Quick Guide
M212843EN	Vaisala Indigo202 Digital Transmitter User Guide
M212816EN	Vaisala Indigo202 Digital Transmitter Quick Guide

1.3 Documentation conventions



WARNING! alerts you to a serious hazard. If you do not read and follow instructions carefully at this point, there is a risk of injury or even death.



CAUTION! warns you of a potential hazard. If you do not read and follow instructions carefully at this point, the product could be damaged or important data could be lost.



Note highlights important information on using the product.



Tip gives information for using the product more efficiently.

1.4 Trademarks

Vaisala® and CARBOCAP® are registered trademarks of Vaisala Oyj.

Modbus® is a registered trademark of Schneider Automation Inc.

Nafion™ is a trademark of The Chemours Company FC, LLC.

Windows® is either a registered trademark or trademark of Microsoft Corporation in the United States and other countries.

All other product or company names that may be mentioned in this publication are trade names, trademarks, or registered trademarks of their respective owners.

2. Product overview

2.1 Introduction to GMP251

GMP251 is designed for CO₂ measurement in demanding applications that require reliable and accurate performance. The measurement range is 0–20 %CO₂.

The probe is based on Vaisala's patented 2nd generation CARBOCAP® technology and equipped with Vaisala Microglow infrared light source. The probe is easy to install with a plug-in/plug-out M12 connection.

Sensor performance is optimized at 5 %CO₂ measurement. For compensation purposes, the probe also includes an internal temperature sensor that allows measurement compensation according to ambient temperature. As dust and most chemicals do not affect the measurement, and the effect of temperature, pressure and background gas can be compensated for, the probe can provide accurate and stable measurements in a wide range of applications.

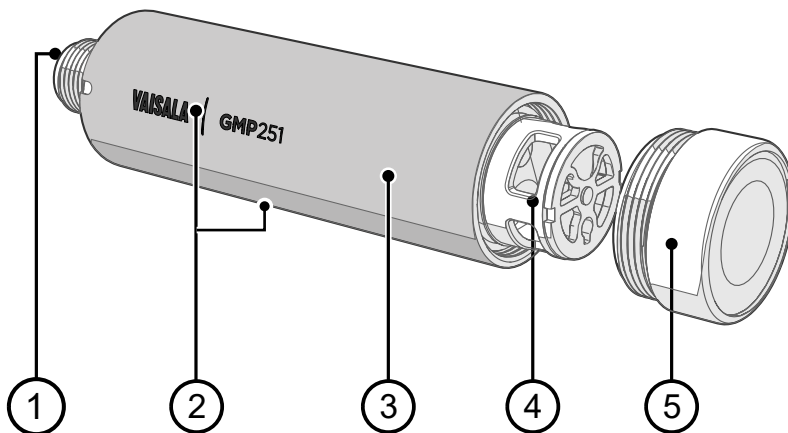


Figure 1 GMP251 probe parts

- 1 5-pin M12 connector. For pinout, see [Wiring \(page 30\)](#).
- 2 Probe name and orientation mark for Vaisala transmitter installations (front) and laser-printed type label (back).
- 3 Probe body. Contains the main component board.
- 4 Measurement cuvette with optics and CARBOCAP® CO₂ sensor.
- 5 Filter (see [Filter options \(page 17\)](#)).



CAUTION! Do not attempt to open the probe body. There are no user serviceable parts inside the probe body.

2.2 Basic features and options

- CO₂ measurement range 0–20 %
- Vaisala CARBOCAP® CO₂ sensor with excellent long-term stability
- Measurement compensated for effects of temperature, pressure, and background gas. The temperature compensation can be based on an integrated temperature sensor or use a set temperature. Pressure and background gas parameters can be set to the probe.
- Heating to avoid condensation on optical elements
- Digital output with RS-485:
 - Modbus RTU
 - Vaisala Industrial Protocol
- Analog output:
 - Current output (0–20 mA or 4–20 mA)
 - Voltage output (0–5 V or 0–10 V)
- Can be used as a standalone probe or with Vaisala Indigo transmitters
- Incorporated into the GMP80P handheld probe, which allows for easy pump-aspirated CO₂ sampling from incubators
- Compatible with Indigo80 handheld indicator and MI70 handheld indicator
- Can be connected to Vaisala Insight PC software for configuration, diagnostics, and temporary online monitoring
- Easy plug-in, plug-out.

More information

- [Operating principle of CO₂ measurement \(page 18\)](#)
- [Environmental compensation \(page 19\)](#)
- [Vaisala Industrial Protocol overview \(page 64\)](#)
- [Modbus \(page 31\)](#)
- [Overview of MI70 support \(page 59\)](#)
- [GMP80P product overview \(page 50\)](#)
- [GMP251 specifications \(page 114\)](#)
- [GMP80P specifications \(page 117\)](#)

2.2.1 Connectivity to Vaisala Insight software

The probe can be connected to Vaisala Insight software using a Vaisala USB cable (item code 242659) or Vaisala Indigo USB adapter (item code USB2). With the Insight software, you can:

- Calibrate and adjust the measurement.
- See device information and status.
- See real-time measurement.

- Configure serial communication settings, purge settings, filtering factor, and analog output parameters and scaling.

More information

- [Connecting to Insight software \(page 32\)](#)

2.2.2 Additional features with Indigo transmitters

GMP251 probes manufactured from 2017 onwards are compatible with Vaisala Indigo200 series transmitters.

Connecting the probe to an Indigo transmitter provides a wide range of additional options for outputs, measurement viewing, status monitoring, and configuration interface access.

Examples of additional features available with Indigo520 transmitters:

- Touchscreen display for real-time data viewing and configuration
- Support for 2 probes simultaneously
- 4 configurable analog outputs
- 2 configurable relays
- 2-wire current loop analog input
- Barometer module for barometric pressure measurement (optional module)
- Ethernet connection with web interface for remote access
- Modbus TCP/IP protocol
- Service port for connecting to Indigo80 handheld indicator

Examples of additional features available with Indigo300 transmitters:

- Numerical and graphical display for real-time data viewing
- 3 pre-configurable analog outputs
- Service port on the front for connecting to Vaisala Insight PC software or Indigo80 handheld indicator for easy configuration access


Examples of additional features available with Indigo200 series transmitters:

- 3.5" TFT LCD color display or non-display model with LED indicator
- Digital output or 3 analog outputs (depending on the transmitter model)
- 2 configurable relays
- USB-C connection to Vaisala Insight PC software for easy configuration access



Available features vary depending on the Indigo transmitter model. For more information on Indigo transmitters, see vaisala.com/indigo.


2.2.3 Filter options

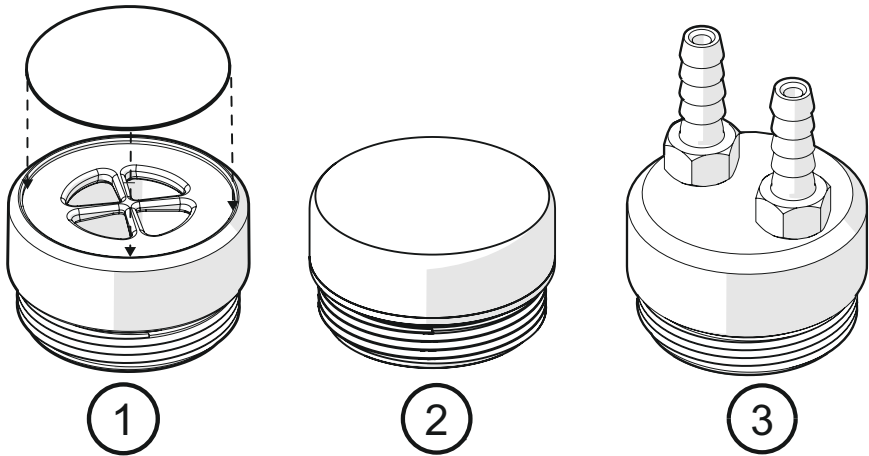

The following filter options are available for GMP251:

1. Standard membrane filter, item code ASM211650SP. Gas can enter only through the top of the filter (plastic grid covered with membrane), the sides of the filter are solid.
2. Porous sintered PTFE filter for extra protection, item code DRW243649SP. The porous material of the filter allows gas to enter from all sides of the filter.
3. Flow-through adapter, item code ASM211697SP. Two gas ports for controlled gas feed (port outer diameter 4.6 mm, port hole inner diameter 2 mm, suitable for tubing with 4 mm inner diameter).

2.3 Operating principle of CO₂ measurement

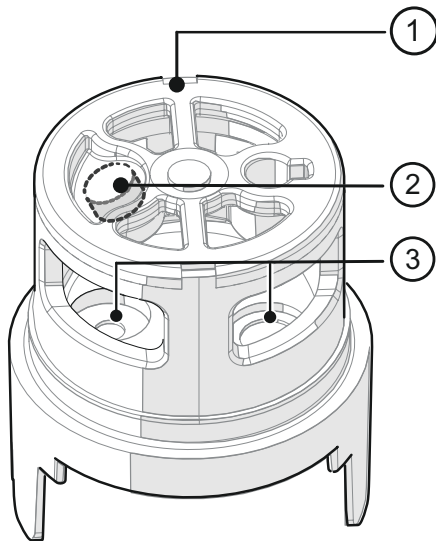

The Vaisala CARBOCAP® sensor used in the probe is a silicon-based, nondispersive infrared (NDIR) sensor for measurement of gaseous carbon dioxide in air-like gases.

Figure 2 Probe cuvette with mirror and sensor chips

- 1 Cuvette
- 2 Mirror
- 3 Sensor chips under TO5 package

The sensitivity to carbon dioxide is based on absorption of infrared light at a characteristic wavelength. During measurement, infrared light is routed through the cuvette that contains the gas to be measured. A mirror reflects the light from the cuvette to a thermopile detector that measures the light intensity at a wavelength determined by a Fabry-Pérot interferometer (FPI) and a band pass filter.

The carbon dioxide measurement consists of two steps: first, the FPI is electrically tuned so that its pass band coincides with the characteristic absorption wavelength of carbon dioxide and the signal is recorded. Second, the pass band is shifted to a wavelength where no absorption occurs in order to get a reference signal. The ratio of these two signals, one at the absorption wavelength and the other at the reference wavelength, gives the fraction of light absorption from which the carbon dioxide concentration is calculated. Measuring the reference signal compensates the possible effects of sensor aging and signal attenuation due to dirt on optical surfaces, making the sensor very stable over time.

TO5 packages with hermetic windows are used to protect the sensor chips from moisture and contamination. A heater chip is utilized to prevent condensation in normal operation.

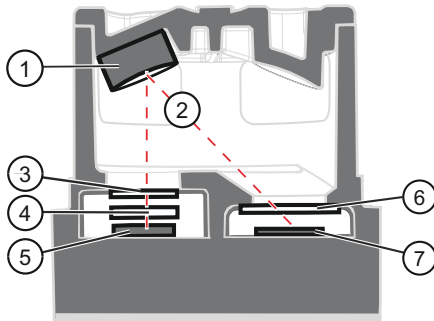


Figure 3 CO₂ measurement in the measurement cuvette

- 1 Gold-plated mirror
- 2 Light absorbed by CO₂ in the measured gas
- 3 Hermetic window
- 4 Fabry-Pérot interferometer
- 5 Light source (Microglow)
- 6 Hermetic window
- 7 Thermopile detector

2.4 Environmental compensation

When necessary, various environmental compensations can be applied to improve the CO₂ measurement accuracy of the probe.

The probe can compensate for the effects of the following parameters:

- Temperature
- Pressure
- Background gas oxygen (O₂) content
- Background gas relative humidity (%RH)



To apply an accurate relative humidity compensation, make sure that also the temperature compensation and pressure compensation configurations match the measurement environment.

The probe has an on-board temperature sensor that can be used to compensate for temperature. Additionally, if the probe is integrated in a system that measures one or more of the compensation parameters (T, P, RH, O₂), they can be updated to the probe continuously.

Compensation parameters are configured on the order form when ordering the probe, and can later be updated using Vaisala Industrial Protocol or Modbus protocol.

You can also turn off any of the compensations. In that case, the probe uses the default compensation value that is mathematically neutral for the probe's internal compensation model.

You can configure the environmental compensation settings with any of the following options:

- Vaisala Industrial Protocol serial commands
- Vaisala Insight software
- Modbus configuration registers
- MI70 handheld indicator
- Indigo80 handheld indicator



The environmental compensation configuration options available in the MI70 handheld indicator are not compatible with the probe. Using MI70 to configure the environmental compensation settings of the probe can decrease measurement accuracy.

More information

- [Environmental compensation commands \(page 88\)](#)
- [Vaisala Insight PC software \(page 32\)](#)
- [Configuring environment settings with Indigo80 handheld indicator \(page 48\)](#)
- [Changing environmental compensation settings with MI70 indicator \(page 62\)](#)
- [Effect of environmental compensations \(page 95\)](#)
- [Configuration registers \(page 127\)](#)

2.4.1 Temperature compensation

The probe can measure the approximate temperature of the CARBOCAP® sensor for compensation, or use a fixed setpoint. Unless a dedicated temperature measurement is available and can be regularly updated to the probe, it is strongly recommended to use the probe's internal temperature compensation to ensure real-time accurate measurements. If the measurement is made in a constant temperature, the constant temperature can be set as the compensation value (fixed setpoint option).

If temperature compensation is turned off, the probe uses the default value of +25 °C (+77 °F).



When the probe is installed through a flange and part of the probe and the cable is left outside the measuring environment, it is possible that temperature conduction from the probe body and cable outside the measurement environment affects the temperature compensation and decreases measurement accuracy.

2.4.2 Pressure compensation

The probe does not have on-board pressure measurement. However, a pressure reading from an external source can be used as a setpoint value for compensation using Vaisala Industrial Protocol or Modbus.

If pressure compensation is turned off, the probe uses the default compensation value of 1013 hPa.

2.4.3 Background gas compensation

The probe does not have on-board oxygen or relative humidity measurement. However, oxygen and relative humidity readings from an external source can be used as setpoint values for compensation via Vaisala Industrial Protocol or Modbus. The default setpoint values are as follows:

- Oxygen concentration: 19.7 %O₂ or 21 %O₂
- Relative humidity: 50 %RH or 93 %RH

If background gas compensations are turned off, the probe uses the value 0 % for both.

More information

- [Environmental compensation \(page 19\)](#)
- [Configuration registers \(page 127\)](#)

2.5 Probe startup

When powered on, the probe starts up within 20 seconds. Measurements from the outputs (digital and analog) become available during this time but note that they will only reach specified accuracy after a 4-minute warm-up period. For this reason, you should design your system so that it does not rely on measurements from the probe during this time.



Specifically note that the CO₂ reading will rise to the correct reading as the sensor's infrared emitter achieves operation temperature.

2.6 Filtering factor

You can set a filtering factor that affects the speed at which the latest CO₂ measurement is integrated into the output of the probe. This allows averaging the output if the measuring environment produces occasional exceptionally high or low readings.

The filtering factor can be set either with the Insight PC software, Modbus configuration registers, or with the Indigo80 handheld indicator.

More information

- [Configuring filtering factor with Insight \(page 41\)](#)
- [Example of configuring filtering factor with Indigo80 handheld indicator \(page 47\)](#)
- [Configuration registers \(page 127\)](#)

2.7 Analog output overrange behavior

The analog output has a defined behavior when the measurement is outside the scaled analog output range. At first, the output is clipped when the measurement exceeds a set limit (the measurement continues, but the output does not change from the clipped value).

When the measurement exceeds the second limit (error limit), the analog output switches to the error state defined for the output.

The table below lists the clipping and error limits and default error state outputs for the analog voltage and current outputs.

Table 3 Analog output overrange clipping and error limits

Output voltage / current	Clipping limit	Error limit	Default error state output
0–5 V	>5 %	>10 %	0 V
0–10 V	>1 %	>10 %	0 V
0–20 mA	>5 %	>10 %	23 mA
4–20 mA	>5 %	>10 %	2 mA

The same clipping and error limits are applied when the measured value drops back to the scaled range: at first the output returns to the clipped value from the error state, and then to normal output.



Clipping and error state limits differ for 0–10 V and 0–5 V outputs. For 0–10 V output the limits are 1 % and 10 %, and for 0–5 V output the limits are 5 % and 10 %.



When setting the error output voltage (**amode** command), note that the probe's maximum voltage output is 10.325 V.

More information

- [Analog output \(page 81\)](#)

2.7.1 Analog output overrange example

Consider a probe with 0–5 V output, scaled to 0–200 000 ppm (= 0–20 %) CO₂.

- When the measured CO₂ rises above 20 %, the output rises above 5 V.
- The output keeps rising until the measurement is 21 %CO₂, at which point the probe outputs 5.25 V.
- If the CO₂ level rises above 21 %CO₂, the output still remains at 5.25 V.
- If the CO₂ level rises above 22 %CO₂, the output enters the error state, which is 0 V for the 0–5 V output.

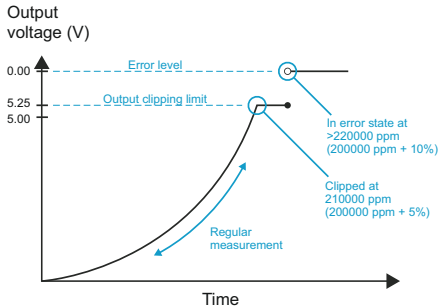


Figure 4 Example of analog output overrange behavior

This example uses output scaled to 0–5 V and 0–200 000 ppm, error level set to 0 V, clipping set to 5 % overrange, and error limit set to 10 % overrange. CO₂ concentrations (ppm) are indicated for the clipping point and error limit point.



This overrange and error behavior is specific to the analog output, and does not affect the readings of the digital outputs.

You can change the analog output overrange behavior using the **aoover** command.

2.8 Safety

The probe delivered to you has been tested for safety and approved as shipped from the factory. Note the following precautions:



WARNING! When returning a product for calibration or repair, make sure it has not been exposed to dangerous contamination, and is safe to handle without special precautions.



CAUTION! Do not modify the unit or use it in ways not described in the documentation. Improper modification or use may lead to safety hazards, equipment damage, failure to perform according to specification, or decreased equipment lifetime.



CAUTION! Do not attempt to open the probe body. There are no user serviceable parts inside the probe body.

2.8.1 ESD protection

Electrostatic Discharge (ESD) can cause immediate or latent damage to electronic circuits. Vaisala products are adequately protected against ESD for their intended use. However, it is possible to damage the product by delivering electrostatic discharges when touching an exposed contact on the product.

To make sure you are not delivering high static voltages yourself, avoid touching the pins on the M12 connector.

3. Installation

3.1 GMP251 probe dimensions

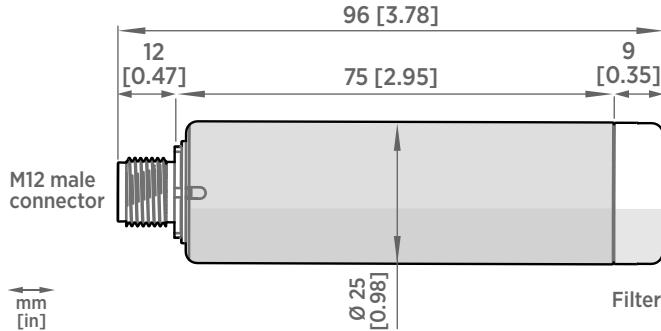


Figure 5 GMP251 dimensions

3.2 Recommended installation

The GMP251 probe can be installed in an environment with an operating temperature range of $-40 \dots +60 \text{ }^{\circ}\text{C}$ ($-40 \dots +140 \text{ }^{\circ}\text{F}$) ¹⁾. Make sure the probe is in a location that represents the measurement environment properly.



In environments with high levels of humidity and risk of condensation, install the probe horizontally to prevent condensation from forming inside the filter and on the optical surfaces.

The 5-pin male M12 connector on the probe provides an easy plug-in/plug-out connection to a compatible cable. Using a shielded cable in the installation and grounding is recommended. The GMP251 probe itself does not require a shielded cable, but if the installation environment contains high levels of EMC noise, the disturbance may be transmitted to another connected device through an unshielded cable. Vaisala Online Store at store.vaisala.com offers a wide selection of shielded cables in different lengths.

More information

- [GMP251 spare parts and accessories \(page 120\)](#)

1) Occasional short-term exposure to up to $+90 \text{ }^{\circ}\text{C}$ ($+194 \text{ }^{\circ}\text{F}$) allowed, provided that the probe is fully installed inside the measured condition and power is switched on. Measurement accuracy specification of the probe not applicable if used in temperatures above $+60 \text{ }^{\circ}\text{C}$ ($+140 \text{ }^{\circ}\text{F}$).

3.3 Installation accessories

The probe can be installed through a surface using the optional flange accessory (Vaisala product code 243261SP), or attached for example to a wall with the optional clip accessory (two-clip set, Vaisala product code 243257SP).

For installations where a pressure washer is used to clean the measurement area, a spray shield with an IPX9 rating (Vaisala product code ASM212017SP) is available as an optional accessory.



In environments with high levels of humidity and risk of condensation, install the probe horizontally to prevent condensation from forming inside the filter and on the optical surfaces.

3.3.1 243261SP mounting flange

The optional flange accessory is used to install Ø 25 mm (0.98 in) probes through a wall or other surface.

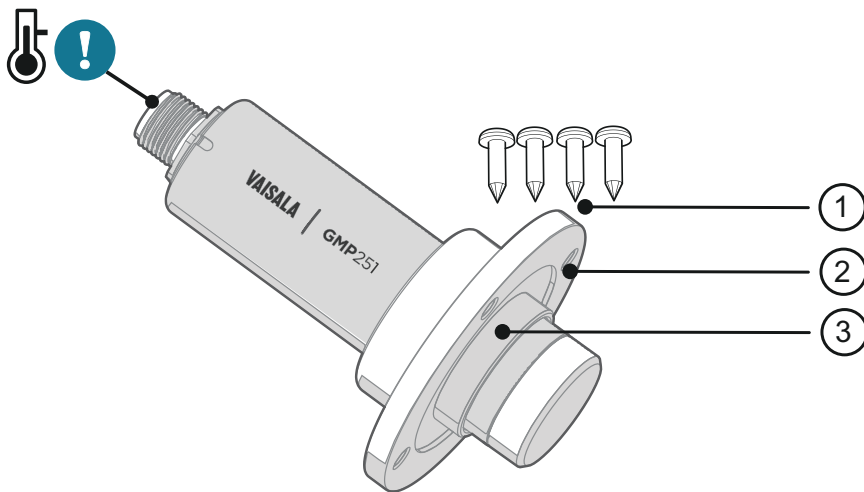


Figure 6 Probe with 243261SP mounting flange

- 1 4 Phillips head screws (included)
- 2 Mounting flange (diameter 60 mm (2.36 in)) with four Ø 4.2 mm (Ø 0.17 in) screw holes
- 3 Gasket ring



Leaving part of the probe body and the cable outside the measurement environment can cause heat conduction that affects the temperature compensation and measurement accuracy.

3.3.2 243257SP mounting clips

The optional mounting clips (set of two clips) are used to hold the probe in place for example on a wall or other surface. Each clip base attaches to the installation surface with one screw (screw hole Ø 4.2 mm).


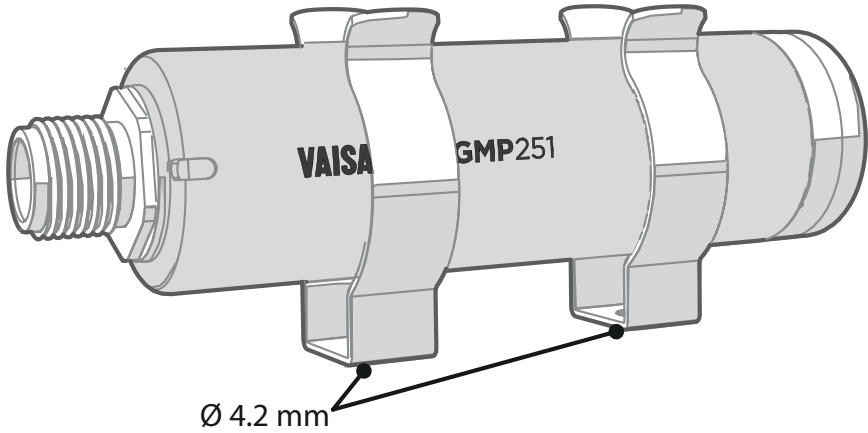


Figure 7 Probe in 243257SP mounting clips

3.3.3 ASM213582 probe holder

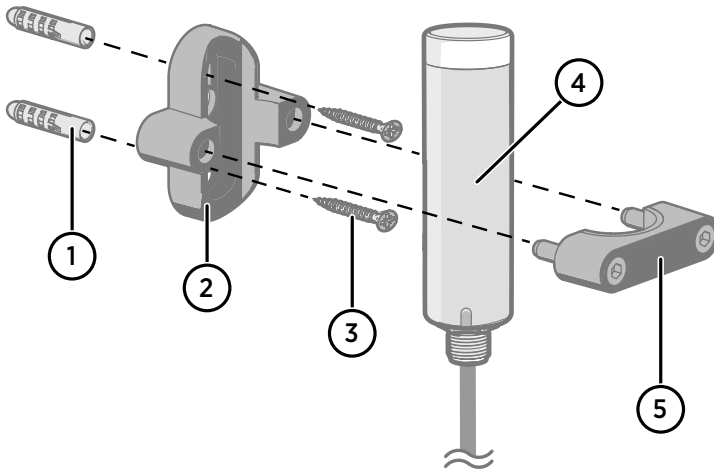

- 1 Wall plugs (2 pcs included, 6×30 nylon)
- 2 Base of the probe holder
- 3 Screws (2 pcs included, 4.8×25 DIN7981C PZ A4)
- 4 Ø 25 mm (0.98 in) probe body
- 5 Top of the probe holder with 2 hex screws (4-mm socket)

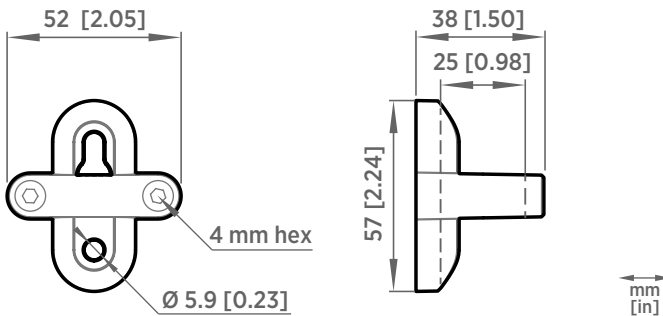


Figure 8 Probe holder ASM213582 dimensions

3.3.3.1 Mounting the probe body using probe holder ASM213582

Probe holder ASM213582 is designed for attaching a 25 mm probe body to a wall or other surface with screws.

- ▶ 1. Use a 4-mm Allen key to open the hex screws and remove the top part of the probe holder.

2. Hold the base part against the mounting surface and mark the locations of the 2 screw holes.
3. Drill holes at the marked locations using a 6-mm drill bit. Make the holes at least 30 mm deep.
4. Insert the wall plugs in the holes.
5. Mount the base part using screws.
6. Insert the probe body in the holder.
7. Place the top part on the base and tighten the hex screws with a 4-mm Allen key.

3.3.4 ASM212017SP spray shield

The optional spray shield (Vaisala item code ASM212017SP) allows washing the probe with a pressure washer without exposing the sensor to moisture (IPX9 rating).

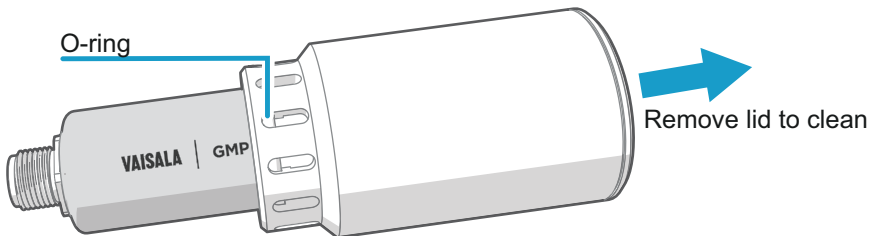


Figure 9 Probe with ASM212017SP spray shield



The spray shield can limit or prevent the use of the probe mounting clips (243257SP).

More information

- [ASM212017SP spray shield dimensions \(page 123\)](#)

3.4 Power supply

The supply voltage range of the probe is 12–30 V DC with the digital output option. If the analog output is used, the supply voltage range is 12–30 V DC for voltage output and 20–30 V DC for current output.

Typical power consumption is less than 0.4 W in continuous operation, and the maximum is 0.5 W.

More information

- [GMP251 specifications \(page 114\)](#)

3.5 Wiring

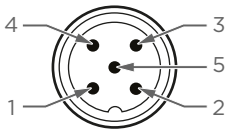


Table 4 M12 male connector

Pin#	Function	Note	Cable 223263SP wire colors
1	Power in	<ul style="list-style-type: none">• With digital output: 12–30 V DC• With voltage output: 12–30 V DC• With current output: 20–30 V DC Typical average power consumption <0.4 W, maximum 0.5 W.	Brown
2	RS-485- or voltage output	Voltage: 0–5 V DC or 0–10 V DC (default analog output scaling)	White
3	GND	–	Blue
4	RS-485 + or current output	Current: 0–20 mA or 4–20 mA (default analog output scaling)	Black
5	Output control	Connecting pin #5 to GND (pin #3) forces the probe to analog output mode. If an analog output configuration has not been selected, default 0–10 V DC and 4–20 mA scalings are used. If pin #5 is not connected, the analog or digital output selected when ordering or set later through configuration is used.	Gray



Note that the probe always remains in analog mode when pin #5 is connected to pin #3, and cannot be switched to digital output in this wiring option.

4. Modbus

The probe can be accessed using the Modbus serial communication protocol. The supported Modbus variant is Modbus RTU (Serial Modbus) over RS-485 interface.

For instructions on enabling the Modbus mode when you are using the probe with Vaisala Industrial Protocol, see [Enabling Modbus mode from Vaisala Industrial Protocol \(page 67\)](#).

For instructions on switching to Vaisala Industrial Protocol when the probe is in Modbus mode, see [Accessing serial commands from Modbus or analog mode \(page 67\)](#).

For information on using Modbus with the Indigo202 digital transmitter, see [Modbus overview \(page 37\)](#).

The pre-configured default Modbus serial settings are presented in the following table.

Table 5 Default Modbus serial communication settings

Description	Default value
Serial bit rate	19200
Parity	N
Number of data bits	8
Number of stop bits	2
Modbus device address	240

More information

- [Modbus registers \(page 126\)](#)

5. Vaisala Insight PC software

Vaisala Insight PC software is a service and utility software for Indigo-compatible Vaisala devices. Insight is available for Microsoft Windows® operating systems (64-bit only).

With the Insight software, you can:

- See device information and status.
- See real-time measurement data.
- Configure serial communication settings, purge settings, filtering factor, and analog output parameters and scaling.
- Calibrate and adjust the device.
- Perform short-term data logging for laboratory testing or evaluation purposes, for example.

Download Insight at vaisala.com/insight.

The probe can be connected to Vaisala Insight software using a Vaisala USB cable (item code 242659) or Vaisala Indigo USB adapter (item code USB2).



GMP251 probes support Insight from probe software version 1.3.0 onwards.

5.1 Connecting to Insight software



- Computer with Vaisala Insight software installed
- USB service cable (item code 242659) or Vaisala Indigo USB adapter (item code USB2)



CAUTION! When connecting several devices at the same time, note that your computer may not be able to supply enough power through its USB ports. Use an externally powered USB hub that can supply >2 W for each port.

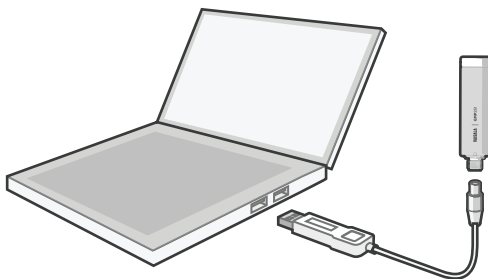


Figure 10 Connecting probe to Insight



Connect the probe to Insight software in the exact order described in the following steps.

- ▶ 1. Open the Insight software on your computer.
- 2. Connect the USB cable to a free USB port on the PC.
- 3. Connect the probe to the USB cable.
- 4. Wait for Insight software to detect the probe.

6. Using GMP251 with Indigo200 series transmitters

6.1 Introduction to Indigo200 transmitters

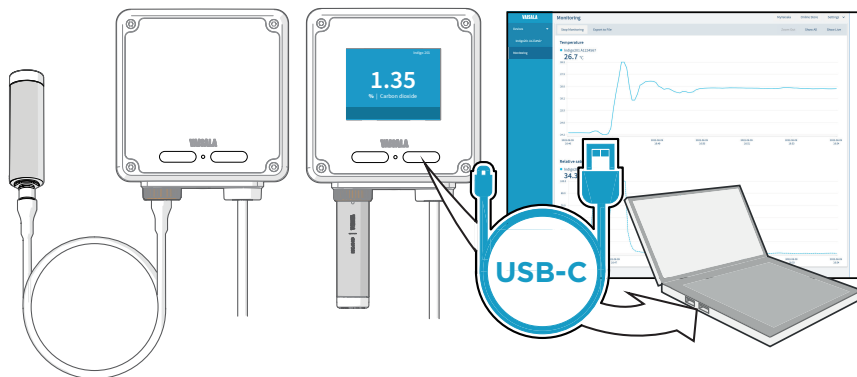


Figure 11 Indigo201 with and without display

Vaisala Indigo200 transmitters are a plug-and-play host device platform for Vaisala Indigo-compatible probes. Indigo transmitters extend the feature set of connected probes with a range of additional options for outputs, measurement viewing, status monitoring, and configuration interface access.

Depending on the Indigo200 transmitter model, a display is available as an optional selection or as a standard feature. In the non-display model, an LED indicator is used for notifications. Probes can be connected either directly on the transmitter from the probe's connector, or by using a cable between Indigo and the probe.

For easy access to configuration and additional monitoring options, the transmitter can be connected to Vaisala Insight PC software with a USB cable.

The Indigo201 transmitter is designed for analog output applications. The transmitter's output options include 3 analog output channels (current or voltage) and 2 relays.

The Indigo202 transmitter is designed for digital output applications. The transmitter's output options include RS-485 Modbus RTU and 2 programmable relays.



Note that Indigo200 transmitters switched from a WLAN configuration interface to a USB-C cable configuration interface (Vaisala Insight PC software) in 2022. For instructions on configuring GMP251 probes with the earlier WLAN interface, see the earlier *GMP251 User Guide* revision M211799EN-G, available at docs.vaisala.com.

For more information on Indigo transmitter models, see vaisala.com/indigo.

6.2 Connecting probes



CAUTION! Only connect compatible measurement devices to the transmitter. Attempting to connect incompatible devices or cables can damage the equipment. Refer to transmitter specifications for compatibility information.

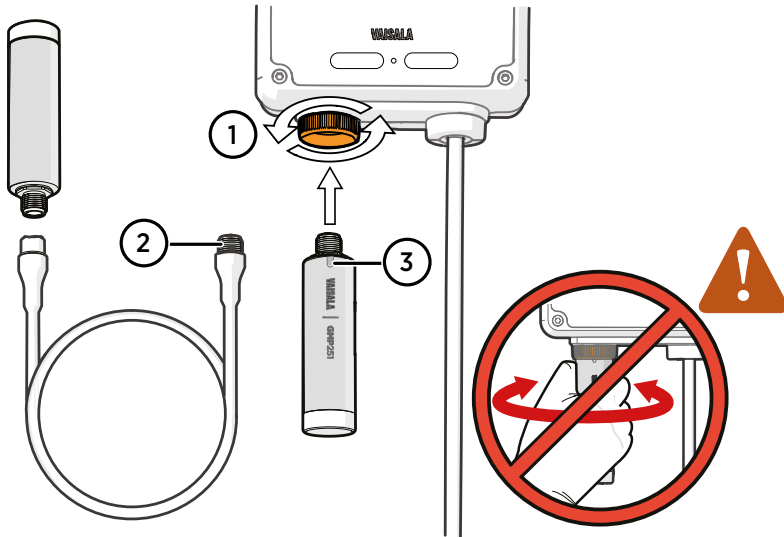


Figure 12 Connecting probes and cables to Indigo200

- 1 Probes are locked in place with the locking wheel. **Never turn from the probe body.**
 - 2 Probe connection cables are connected in the same way as probes: insert the cable in the connector and hold in place while turning the locking wheel.
 - 3 Probes are inserted into the probe connector with the orientation mark and/or the Vaisala logo facing out.
- 1. Insert the probe into the probe connector with the orientation mark and/or the Vaisala logo on the probe body facing out. Note the orientation of the cable connector pins when inserting.
2. Hold the probe in the probe connector and lock it in place by turning the locking wheel **counterclockwise**. Never turn the probe body when connecting or removing probes, only the locking wheel of the transmitter.

- 3. If no previous configuration exists (that is, this is the first time a probe is connected to the transmitter, and no configuration has been set manually), the transmitter adapts the analog output settings of the connected Vaisala Indigo-compatible probe automatically.



If a configuration that differs from the connected probe's configuration exists in the transmitter, you must configure the new probe manually to enable the analog output.

- 4. When the transmitter recognizes the connected probe, it shows a notification message on the display (for example, **Probe Connected: GMP251**). In the non-display model, a green LED is lit when the transmitter recognizes the probe.

6.3 Analog output configuration overview

The transmitter has 3 analog outputs that can be configured with the Insight PC software. Each output has the same set of configuration options in the analog output 1–3 menus under the **Configure device** menu (**Analog output 1** shown in example).



The available selections depend on the probe that is connected to the transmitter (for example, CO₂ output parameter selection is only available for CO₂ probes).

Configure Device

Indigo201 A1234567

In Advanced Mode

General

Display

Displayed parameters

Displayed units

Analog output 1

Analog output 2

Analog output 3

Relay output A

Output mode

Output parameter

Scale low end

Scale high end

Error output level

Clipping limit

Error limit

4...20 mA

CO₂

100

500

3.60

1.0

2.0

mA

%

%

Save

Close

Figure 13 Analog output configuration options example in Insight

The **Analog output type** selection in the **General** menu defines whether the outputs are used for voltage or current output. When you change the output type, Insight reloads the transmitter.

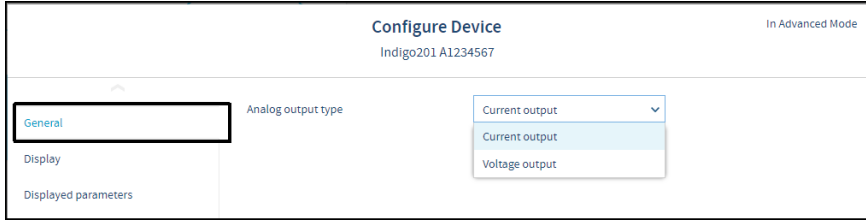


Figure 14 Analog output type (current or voltage) selection in Insight

Receiving analog output settings automatically from probe

The transmitter can receive the analog output settings automatically from the connected probe in the following cases:


1. The transmitter is a new unit and has not been configured yet.
2. The user has removed the analog output configuration from the transmitter.



It is recommended to always review the analog output settings with the Insight PC software after connecting a probe.

6.3.1 Analog output configuration example

- 4–20 mA current output

1. Connect the transmitter to the Insight software.
2. Select  and open the **Configure device** menu.
3. Open the **General** menu and select **Current output** from the dropdown. Store your selection with **Save**.
4. Return to the **Configure device** menu and select one of the 3 analog outputs.
5. Select **4 ... 20 mA** from the **Output mode** dropdown.
6. Type **0** in the **Scale low end** field.
7. Type **3.6** in the **Error output level** field.
8. Select **Save** to take your configuration into use.

6.4 Modbus overview

Indigo202 supports the Modbus RTU serial communication protocol (2-wire RS-485 interface).

There are 2 groups of Modbus register addresses in use in Indigo202: probe registers and Indigo registers. The probe registers are received from the connected probe, and are organized according to the probe's register map. Indigo registers include transmitter-specific information.

Table 6 Indigo202 Modbus registers

Address	Name	Data type
Probe registers (according to the connected probe's register map) ¹⁾		
0x0000	First measurement probe address	
0xDFFF	Last measurement probe address	
Indigo registers		
0xE000	Status	16-bit
0xE001	Notification and error bits	16-bit
0xE002	Connected probe	text [30]
0xE011	Relay A status	enum
0xE012	Relay B status	enum

1) See the connected probe's Modbus documentation for probe-specific register information

Received Modbus requests for register operations are handled in two different ways, depending on whether the register is a probe register or an Indigo transmitter register. Addresses above 0xDFFF (that is, Indigo registers) are handled as normal requests. Lower addresses (that is, probe registers 0x0000–0xDFFF) are passed to the measurement probe, and the response from the probe is again passed to the original Modbus client. Indigo can also have a cache for commonly requested registers (Measurement registers).

The maximum response delay is 2 seconds (when the data content needs to be fetched from the probe). The minimum delay between requests is 10 ms.

6.4.1 Default communication settings

Figure 15 Insight PC software Communication menu view

Modbus communication settings can be configured in the Insight PC software **Configure device > Communication** menu.

Table 7 Default Modbus serial communication settings

Description	Default value
Serial bit rate	19200
Parity	None
Number of data bits	8
Number of stop bits	2
Modbus device address	10

6.5 Relay configuration overview

Indigo200 transmitters have 2 configurable relays (relay A and relay B). Both relays have configuration options for selecting the parameter that is used to control the relay, activation triggers, hysteresis, and error state behavior.

Relay output on/off	<input checked="" type="checkbox"/>	1
Output parameter	CO ₂	2
Activate above	500	3
Activate below	100	4
Hysteresis	20	5
Error output state	Inactive	6

Figure 16 Relay configuration options (CO₂ probe example)

- 1 **Relay output on/off:** Set to **ON** to enable relay output.
- 2 **Output parameter:** The measurement that is used to control the relay. Options vary based on the type of the connected probe.
- 3 **Activate above:**
When the measurement value exceeds the value set here, the relay activates.
- 4 **Activate below:**
When the measurement value falls below the value set here, the relay activates.
- 5 **Hysteresis:** Define a buffer value for relay deactivation (set the value according to the parameter selected for the relay). With hysteresis, an activated relay switches off only when the measured value crosses the relay activation/deactivation limit and the additional buffer value.
- 6 **Error output state:** Select which state the relay is set to when an error occurs (on, off, or remains in its current state).

Relay state information

When one or both relays are enabled, the relay state (active/not active) is shown on the transmitter display. You can also view the status of the relays with Insight PC software.

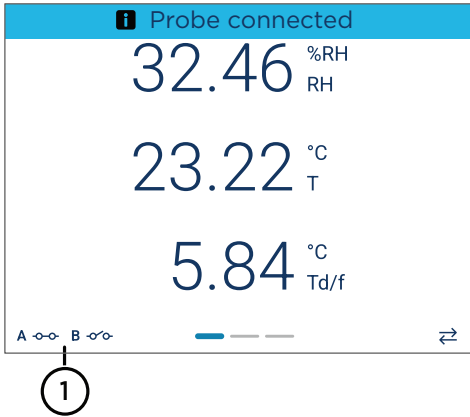


Figure 17 Relay icons on display (relay A active, relay B not active)

6.6 Configuring filtering factor with Insight

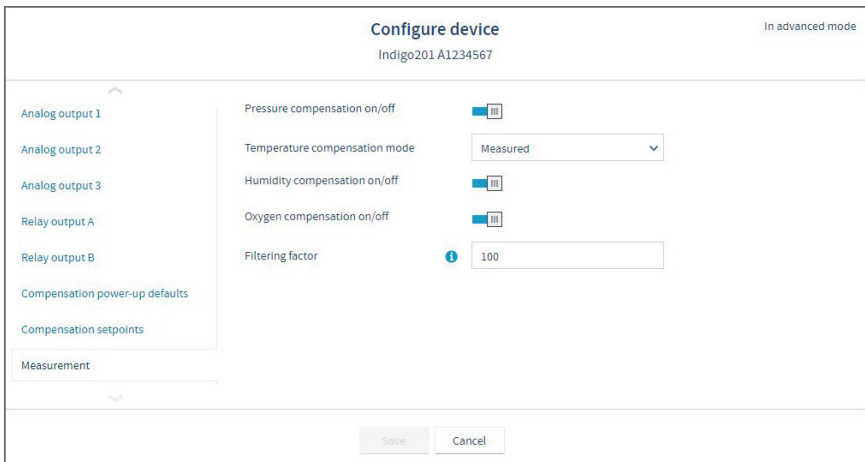



Figure 18 Filtering factor setting in Insight

The filtering factor affects the speed at which the latest measurement is integrated into the output of the probe: a new measurement is produced approximately every two seconds. You can configure the filtering factor in Insight PC software by selecting  > **Configure device** > **Measurement**. The option is available only in **Advanced mode** of Insight.

By default, the filtering factor is set to 1.0 (1.0 = 100 in the configuration range of the Insight user interface), which means the latest measurement is shown directly in the output, without any filtering. To apply filtering, enter a lower filtering factor to include previous measurements in the calculation of measurement output. For example, changing the filtering factor to 0.1 results in an output that is a combination of previous measurements (90 %) and the latest measurement (10 %).

The configuration range of the filtering factor is 0–100: for example, to set the factor to 0.5, set the value to 50 in the **Filtering factor** field.

More information

- [Filtering factor \(page 21\)](#)

7. Using GMP251 with Indigo80 handheld indicator

7.1 Indigo80 handheld indicator

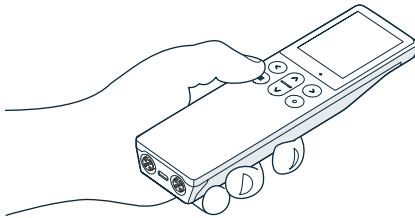


Figure 19 Indigo80 handheld indicator


Vaisala Indigo80 Handheld Indicator is a portable diagnostics tool that accommodates up to two Vaisala Indigo-compatible probes or transmitters for measuring a wide range of parameters.

With the indicator, you can:

- See real-time measurements and device and status information
- Log measurement data
- Calibrate and adjust the probe
- Configure probe features and settings such as compensation setpoints, filtering factor, and serial communication. The available features and settings depend on the probe model and firmware version.



Accessing certain configuration options for your probe requires using the free Insight PC software, downloadable at vaisala.com/insight.

The help tours in the indicator's user interface guide you through the key features of the indicator. You can access the tours in the **Help** menu by pressing the  button.

For more information on using the indicator, for example, editing the measurement views and performing data logging, see [Indigo80 User Guide \(M212722EN\)](#).

More information

- [Configuring probe features with Indigo80 handheld indicator \(page 47\)](#)
- [Configuring environment settings with Indigo80 handheld indicator \(page 48\)](#)

7.1.1 Probe compatibility

The Indigo80 handheld indicator is tested for compatibility with probes that have firmware version 1.3.6 or newer.

Probes with older firmware versions may have limited compatibility with the indicator.

For the most up-to-date version compatibility information, see [Indigo80-compatible Firmware Technical Note \(M212901EN\)](#).

7.1.2 Indigo80 keypad

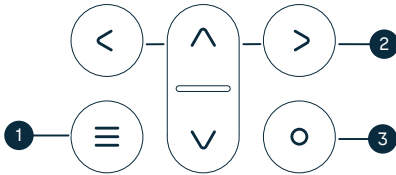



Figure 20 Indigo80 keypad

- 1 Power on/off and main menu button
- 2 Arrow buttons for navigating menus and scrolling views
- 3 Select button for selecting items in the user interface

7.1.3 Indigo80 main menu

Pressing the  button while navigating the Indigo80 menus and views opens the main menu.

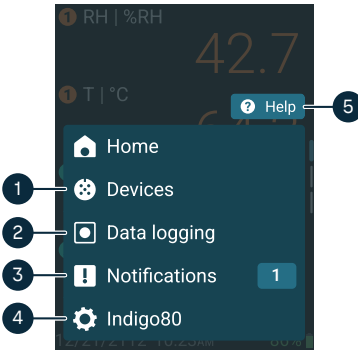


Figure 21 Indigo80 main menu

- 1 **Devices** menu contains, for example, options related to sensor purge, calibration, and environment settings (depending on the connected device).
- 2 **Data logging** menu for setting logging interval and duration, and browsing data files.
- 3 **Notifications** menu displays notifications related to Indigo80 and the connected devices.
- 4 **Indigo80** menu for changing the settings of Indigo80 (for example, date, time, and language) and viewing device information.
- 5 **Help** menu contains tours showing the key features of Indigo80, as well as instructions for sending devices to Vaisala for calibration and maintenance.

7.2 Connecting probes to Indigo80



- Probe connection cable

Up to two Vaisala Indigo-compatible probes or transmitters can be connected to the ports located on the bottom of Indigo80. You can connect and disconnect devices both when the indicator is powered on and when it is off.

Vaisala recommends using cables provided by Vaisala when connecting devices to the indicator. Cables and other accessories are available to order at store.vaisala.com.

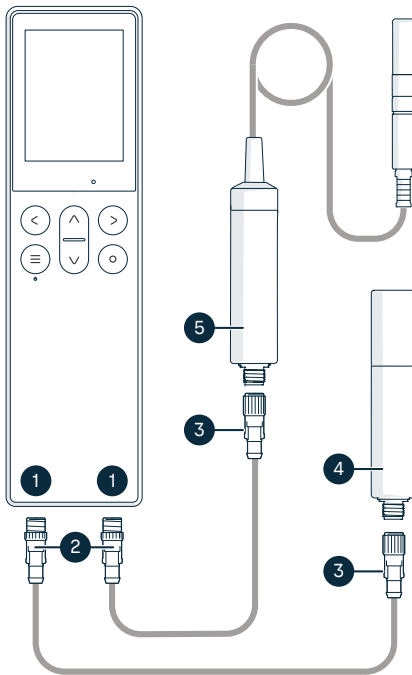



Figure 22 Connecting probes to Indigo80

- 1 M12-5F ports on the bottom of Indigo80 for connecting compatible Vaisala devices. Ports are labeled ① and ②.
- 2 M12-5M A-coded cable connector
- 3 M12-5F A-coded cable connector
- 4 Probe displayed as ① by Indigo80 (GMP252 shown)
- 5 Probe displayed as ② by Indigo80 (HMP7 shown)


- 1. If the indicator is powered and no devices are connected to it, the text **Please connect a probe** will be shown on the display.

- Insert the probe connection cable in one of the ports on the bottom of the indicator.

- **Note the orientation of the cable connector when inserting it**
- **Hold the connector in place while turning its locking ring clockwise – never twist the connector body!**

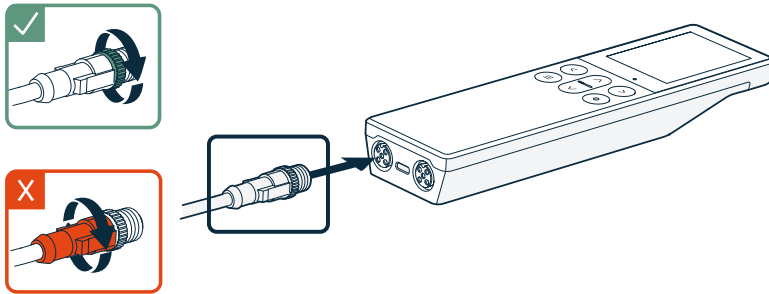

- Connect the probe to the M12-5F end of the probe connection cable.

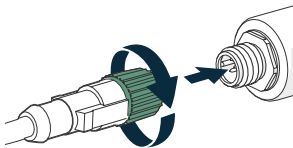

When the indicator recognizes the connected probe, it shows a notification on the display (for example, **GMP252 connected**). A probe connected to the leftmost port in the indicator is labeled ❶ on the indicator's display, while the probe in the rightmost port is labeled ❷.

- To change probes, simply detach the cable from the probe and connect a new probe.



For optimal measurement accuracy, the indicator guides you to check the connected probe's environment settings next.

More information

- [Connecting GMP80P to Indigo80 \(page 51\)](#)

7.3 Configuring probe features with Indigo80 handheld indicator



- Probe connection cable

You can use the Indigo80 handheld indicator to configure probe features such as compensation setpoints, filtering factor, and serial communication. The available features depend on the probe model and firmware version.



Accessing certain configuration options for your probe requires using the free Insight PC software, downloadable at vaisala.com/insight.

1. Connect the probe(s) to the indicator.
2. Open the indicator main menu by pressing
3. Select **Devices**. If you have more than one device connected to the indicator, make a further selection between the devices.
4. Select **Settings** to access and change the features available for your probe.
5. Exit the menu by pressing or return to the main menu by pressing

More information

- [Environmental compensation \(page 19\)](#)
- [Filtering factor \(page 21\)](#)

7.3.1 Example of configuring filtering factor with Indigo80 handheld indicator






- Probe connection cable

The filtering factor allows averaging the output if the measuring environment produces occasional exceptionally high or low readings.

This example shows how to configure the filtering factor as **90** using the Indigo80 handheld indicator.

This means that the new measurement output is a combination of the latest measurement (90 %) and the previous measurement (10 %).

1. Connect the probe(s) to the indicator.

2. Open the indicator main menu by pressing .
3. Select **Devices**. If you have more than one device connected to the indicator, make a further selection between the devices.
4. Select **Settings > Measurement**.
5. Select **Filtering factor**.
 - a. Use the arrow buttons to set the value as **90**.
 - b. Select **Set**.
6. Exit the menu by pressing  or return to the main menu by pressing .

More information

- [Filtering factor \(page 21\)](#)

7.4 Configuring environment settings with Indigo80 handheld indicator



- Probe connection cable


For optimal measurement accuracy, you can use the Indigo80 handheld indicator to give additional information about the connected probe's measurement environment.

For example, you can:

- Configure the pressure compensation setting for your probe
- Enable temperature compensation from another Indigo-compatible probe that has a temperature output



The environment settings available through Indigo80 correspond to the environmental compensation settings of the probe. Therefore, **Environment settings** and **Settings > Compensation setpoints** are interconnected: the configuration set in either menu is applied to both.

- ▶ 1. Connect the probe(s) to the indicator.
- 2. Open the indicator main menu by pressing .
- 3. Select **Devices**. If you have more than one device connected to the indicator, make a further selection between the devices.

4. Select **Environment settings** to access the settings for your probe.

The available environment settings depend on the probe and probe firmware version.

Pressure

Temporary value for pressure compensation. At device restart, compensation value changes back to power-up default value.

Temperature

Temporary value for temperature compensation. At device reset, compensation value changes back to the power-up default value.

Relative humidity

Temporary value for relative humidity compensation. At device reset, compensation value changes back to the power-up default value.



Oxygen concentration

Temporary value for oxygen concentration compensation. At device reset, compensation value changes back to the power-up default value.

5. Select a setting from the list, for example, **Pressure**.

- a. Select the unit, if prompted.

- b. Use the arrow buttons to set the desired value, and then select **OK** or **Set**.

You can also select to read the value from another probe connected to the indicator, provided that the probe can output the required compensation parameter. Press  to highlight the **Read value from** field, and then press  to select the other probe as the source of the value.

6. Scroll down the **Environment settings** view and check that the compensation settings corresponding to the environment settings you defined in the previous steps are set to **On**.

For example:

- If you set a value for **Pressure**, make sure that **Pressure compensation on/off** is set to **On** as well.
- If you set a value for temperature, make sure that **Temperature** is set to **Setpoint**.

7. Exit the menu by pressing  or return to the main menu by pressing .

8. GMP80P Portable CO₂ Probe with Pump Sampling

8.1 GMP80P product overview



Figure 23 Vaisala GMP80P CO₂ probe

Vaisala GMP80P CO₂ probe, which incorporates the GMP251 probe, has been designed for portable sampling and spot-checking use with incubators and fixed CO₂ transmitters.

CO₂ measurement data can be monitored with Vaisala Indigo80 Handheld Indicator connected to GMP80P. Indigo80 can be used for short-time logging of measurement data, as well as for calibrating and adjusting GMP80P. GMP80P is powered by Indigo80, which can supply power for up to two GMP80P probes.

GMP80P can also be connected to a PC running Vaisala Insight PC software for Windows®, which can be used for field calibration, device analytics, configuration functionality, and powering of GMP80P.

More information

- [Vaisala Insight PC software \(page 32\)](#)
- [Indigo80 handheld indicator \(page 43\)](#)
- [GMP80P specifications \(page 117\)](#)
- [GMP80P dimensions \(page 122\)](#)

8.2 GMP80P probe structure

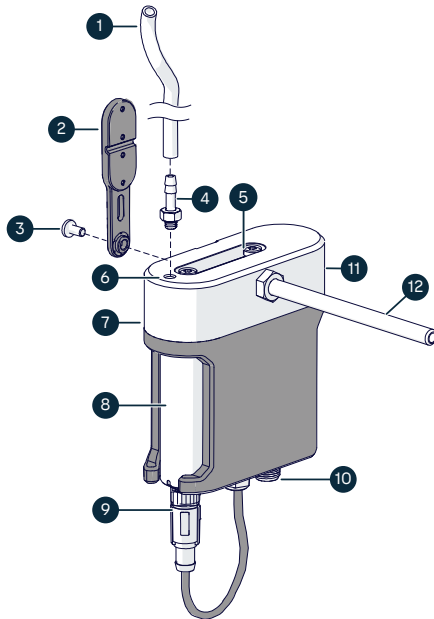


Figure 24 GMP80P probe structure

- 1 Optional plastic tube for gas outlet (not provided by Vaisala)
- 2 Magnetic hanger for attaching GMP80P to metallic surfaces
- 3 Hex-head screw (M4×8) for magnetic hanger
- 4 Gas port, suitable for ID 4 mm (approx. ID 5/32 in) polyurethane or silicone tubing
- 5 Condensation window
- 6 Gas outlet
- 7 CO₂ measurement chamber
- 8 Vaisala GMP251 CO₂ probe
- 9 M12 cable connector for CO₂ probe
- 10 M12 connector for probe connection cable 272075SP (connection to Indigo80 indicator)
- 11 Type label at the back of the device
- 12 Gas inlet through steel pipe (optional polyurethane or Nafion™ membrane tube also available)

8.3 Connecting GMP80P to Indigo80



- Probe connection cable (Vaisala item code 272075SP)



For instructions on connecting GMP80P to the Indigo80 indicator, see [Connecting probes to Indigo80](#) (page 45).


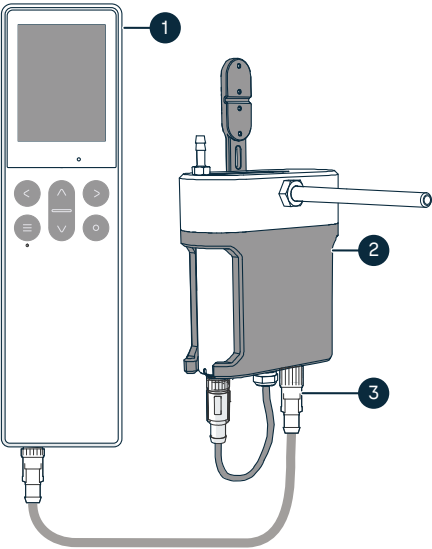


Figure 25 GMP80P connected to Indigo80

- 1 Indigo80 indicator
- 2 GMP80P probe
- 3 Probe connection cable 272075SP connected to the rightmost M12 connector of GMP80P

Note that Indigo80 identifies the GMP80P probe as **GMP251**:



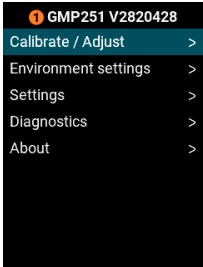



More information

- [Connecting probes to Indigo80 \(page 45\)](#)

8.3.1 Configuring GMP80P for sampling use

Before starting measurements, configure the GMP251 probe in GMP80P for sampling use with Indigo80:

- ▶ 1. **Set temperature compensation mode:**¹⁾
 - a. With GMP80P connected to Indigo80, open the indicator main menu by pressing .
 - b. Select **Devices**. If you have more than one device connected to the indicator, make a further selection between the devices.
 - c. Select **Settings**.
 - d. Select **Measurement > Temperature compensation mode**.
 - e. Select **Measured** (if not already selected).
 - f. Press  to return to the **Settings** menu.
2. **Set relative humidity compensation value:**
 - a. In the **Settings** menu, select **Compensation setpoints > Relative humidity**.
 - b. Depending on your sampling method, enter the relative humidity value as follows:
 - **Steel pipe:** 97 %RH
 - **Plastic tube:** 90 %RH
 - **Nafion tube combined with plastic tube:** relative humidity value of ambient air in the location where your GMP80P is located (for example, 40 %RH)
 - c. Select **Set** to take the value into use.
 - d. Exit the menu by pressing  or return to the main menu by pressing .

8.4 GMP80P installation and sampling

Typical conditions inside incubators are warm and humid (37 °C (98.6 °F), > 90 %RH). This means that condensed water will form easily in the tubing between GMP80P and the incubator, or inside the measurement chamber of GMP80P (when sampling with the steel pipe). To prevent condensation from affecting measurement accuracy, note the condensation prevention considerations in sections [Sampling using plastic tube \(page 57\)](#) and [Sampling using steel pipe \(page 58\)](#).

Before starting measurements, see section [CO₂ dilution coefficients for pump sampling method \(page 54\)](#). The section contains tables to help you determine the most accurate measurement values in your sampling environment.

8.4.1 Nafion membrane tube in sampling

For condensation prevention in sampling, Vaisala recommends using a Nafion membrane tube to maximize transfer of water vapor from the sample to ambient air. The Nafion material is highly selective in the removal of water. Water moves through the membrane wall and evaporates into the surrounding air in a process called pervaporation. For more information about the technology, see permapure.com.

1) For GMP251 probes with configuration code A1C0C0N1, the correct temperature compensation mode has already been set at the factory.



Figure 26 Nafion tube, Vaisala item code 212807SP

The following accessories for tube sampling are available, either when ordering the GMP80P probe or from store.vaisala.com:

Description	Item code
Sampling tube (Nafion membrane and polyurethane), 2.1 m (6 ft 11 in) ¹⁾	GMTUBESP
Sampling tube (Nafion membrane), 0.6 m (1 ft 12 in)	212807SP
Sampling tube (polyurethane), 1.5 m (4 ft 11 in)	279471SP

1) Consists of a 0.6-m (1 ft 12 in) Nafion tube and a 1.5-m (4 ft 11 in) polyurethane tube.



To prevent breaking the Nafion tube, do not insert it directly through the incubator sampling port or attach it to GMP80P. Instead, use a short section of plastic tube (for example, polyurethane tube) at both ends of the Nafion tube. For a setup example, see [Figure 27 \(page 57\)](#).

Plastic tube specifications:

- Inner diameter: 4 mm (approx. 5/32 in)
- Outer diameter: 6 mm (approx. 15/64 in)

Vaisala spare part item GMTUBESP consists of a Nafion tube and a compatible polyurethane tube for this purpose.

8.4.2 CO₂ dilution coefficients for pump sampling method

When using drying tubing (for example, Nafion) for sampling, the CO₂ concentration of the dry sample will be slightly higher than in the wet sample. This is due to water condensing off the warm, humid gas sample as it cools down. In a sampling environment with high dew point temperature, the CO₂ concentration is diluted by the volume that the water vapor occupies. If water vapor is removed from the sample, the fractions occupied by other gases, including CO₂, will increase accordingly.

Table 8 (page 55) contains the dilution coefficients for the gas concentration when drying a gas sample using Nafion tubing. The ranges shown in the table represent a typical environment for GMP80P measurements. Table 71 (page 134) in the appendix of this document shows the data with a wider dew point temperature scale.

When you have the CO₂ measurement value shown by Indigo80, you can use the coefficients given in the tables to calculate the actual gas concentration in the sampling environment:

- ▶ 1. On the horizontal axis of the table, select the dew point temperature of the sampling environment (for example, an incubator).
2. On the vertical axis, select the dew point temperature of the environment where GMP80P is located.
3. Using the coefficient given for these temperature values in the table, multiply the CO₂ value shown by Indigo80.

For example:

- T_d of sampling environment: **+40 °C T_d**
- T_d of environment where GMP80P is located: **+10 °C T_d**
- CO₂ concentration shown by Indigo80: **5.32 %CO₂**
- Actual CO₂ concentration at sampling environment: **5 %CO₂**
(5.32 %CO₂ × 0.939 = 5 %CO₂)

Table 8 CO₂ dilution coefficients for measurements with Nafion tube, T_d of sampled gas 0 ... +60 °C

	Dew point temperature of sampled gas (°C)						
T _d (°C) ¹⁾	0	+10	+20	+30	+40	+50	+60
-20	0.995	0.999	0.978	0.959	0.928	0.879	0.804
-10	0.997	0.990	0.979	0.961	0.930	0.881	0.806
0	1.000	0.994	0.983	0.964	0.933	0.884	0.809
10		1.000	0.989	0.970	0.939	0.890	0.815
20			1.000	0.981	0.950	0.901	0.826

1) Dew point temperature of environment where GMP80P is located.



You can use the Vaisala Humidity Calculator (available at [vaisala.com](https://www.vaisala.com)) to determine the dew point temperature of your sampling environment when you know its temperature and relative humidity values.

Table 9 (page 56) lists the three sampling methods of GMP80P and gives the dilution coefficients for each type of method.

Table 9 CO₂ dilution coefficients for GMP80P sampling methods at 25 °C (77 °F) ambient temperature

	Dew point temperature of sampled gas (°C)								
Method	+20	+25	+30	+35	+40	+45	+50	+55	+60
Nafion + plastic tube (T _d =10 °C)	0.989	0.980	0.970	0.955	0.939	0.915	0.890	0.853	0.815
Plastic tube	0.999	0.995	0.988	0.978	0.965	0.949	0.930	0.908	0.883
Steel pipe	1.000	1.000	0.997	0.991	0.982	0.970	0.955	0.937	0.916

8.4.3 Sampling using plastic tube

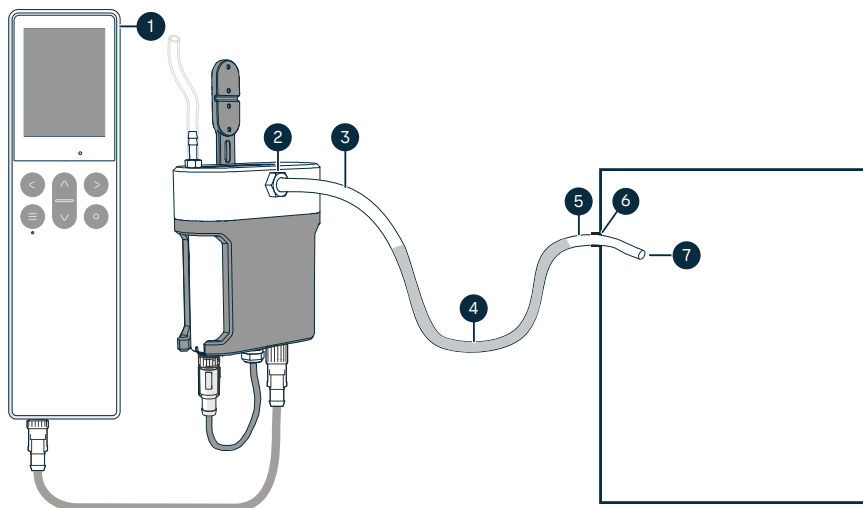


Figure 27 CO₂ sampling from incubator using tube, with GMP80P connected to Indigo80

- 1 Indigo80 indicator
- 2 Adapter for sampling tube (Vaisala item code 279473SP)
- 3 Sampling tube (polyurethane or silicone), OD 6 mm (approx. 15/64 in)
- 4 Section of Nafion membrane tube (moisture exchanger) for optimal condensation prevention. Recommended minimum length approx. 25 cm (10 in).
- 5 Short section of silicone or polyurethane tube. Nafion membrane tube is not recommended to be inserted through the incubator sampling port.
- 6 Incubator sampling port
- 7 Tube inserted into the measurement environment by approx. 5 cm (2 in)

- Position GMP80P higher than the incubator sampling port and leave a loose length of tube hanging between the probe and the incubator (see the figure above).
- Make sure the tube is inserted fully (by approx. 5 cm (2 in)) into the measurement environment.
- Vaisala recommends using a Nafion membrane tube for sampling to maximize transfer of water vapor from the sample to ambient air. For more information about the tube, see [Nafion membrane tube in sampling \(page 53\)](#).
- There is no specific time limit for sampling, but **you need to remove condensed water from the tube from time to time**. Condensation will start to form inside the tube within 5 minutes, actual droplets of water within hours.

To remove condensed water from the sampling tube:

- ▶ 1. Stop the measurement by switching off Indigo80 or by disconnecting GMP80P from Indigo80.

2. Detach the tube from GMP80P and the incubator.
3. Shake the tube to remove any condensed water from inside. Alternatively, you can use compressed air to blow moisture out from the tube.
4. Leave the tube to dry in room temperature.
5. Reattach the tube to GMP80P and the incubator sampling port and continue the measurement.

8.4.4 Sampling using steel pipe

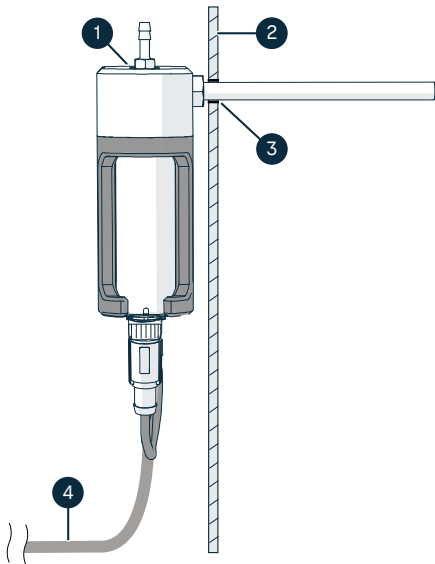


Figure 28 CO₂ sampling from incubator using pipe

- 1 Condensation window
- 2 Incubator wall
- 3 GMP80P sampling pipe inserted through the incubator sampling port
- 4 Connection cable to Indigo80

- Recommended maximum sampling time with the pipe is 15 minutes. Monitor the level of condensation through the condensation window on top of GMP80P.
- Let GMP80P dry between measurements. For example, after a sampling time of 5 minutes, 10 minutes of drying time is required.

To remove condensed water from GMP80P:

1. Detach GMP80P from the incubator wall (do not disconnect from Indigo80).
2. Place GMP80P on a flat surface and leave it running, making sure it does not fall off the surface.

You can continue measuring when there is no more condensation visible in the window on top of GMP80P.

9. Operating with MI70 indicator

9.1 Overview of MI70 support

The GMP251 probe is compatible with instruments that utilize the MI70 indicator, for example the GM70 Handheld Carbon Dioxide Meter. The MI70 indicator is a convenient service tool for viewing the measurement readings and performing calibration and 1-point adjustment.

9.1.1 Probe warm-up and MI70 automatic power off

Note that the probe's warm-up time to full accuracy specification when starting up the probe can be several minutes. When you power the probe with the MI70 indicator, disable the automatic power off feature of MI70 to keep the probe powered on continuously in order to avoid probe restart. To ensure powering during measurements, it is recommended to use the MI70 charging adapter when possible.

9.2 Basic display

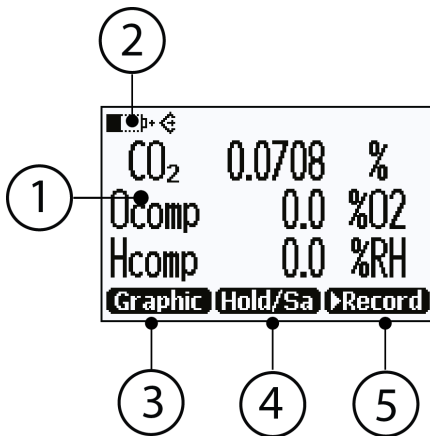


Figure 29 MI70 basic display

- 1 Measured parameter and compensations (up to three items on display simultaneously). You can change the shown items in **Main menu > Display > Quantities and units**.
- 2 Battery indicator. Shows current status (charge) of the battery.
- 3 Function key **Graphic** shows the readings as a curve.
- 4 Function key **Hold/Save** freezes the display and you can save the reading in the MI70 memory.
- 5 Function key **Record** is a quick access to the **Recording/Viewing** menu.

You can change the default function key shortcuts (**Graphic**, **Hold/Save**, **Record**) to other menus or functions in **Main menu > Settings > User interface > Program shortcut keys**.

9.3 Graphical display

The graphical display shows you the measurements as a curve (the curve of the uppermost parameter shown in the basic display). From the curve you can examine the data trend and history of the last minutes.

To open the graphical display, select **Graphic** in the basic display or select **Main menu > Display > Graphic history > Show**.

To get the statistical info on the graph area (minimum, maximum, and average values), press **Info**.

To get the curve of the other selected parameters, press **Next**. To get the curves of all the parameters, press **Next** until the text **All** appears, and then select **All**.







To zoom in and out, press the up/down arrow buttons.

To move back and forward in the timeline, use the left/right arrow buttons.

9.4 Main menu

In the main menu, you can configure the MI70 settings and basic display options, view information about the probe, access recordings and clear the memory, set alarms, start adjustments, and use the analog output option of the MI70 indicator.

To open the main menu and navigate in the menus:

- ▶ 1. Go to the basic display.
2. Press any arrow button, then select  **Open** (must be pressed within 5 seconds or the indicator returns to the basic display).
3. Move in the menu using the   buttons.
4. Select an item with the  button.
5. To return to the previous level, press .
6. To return to normal operation, press  **Exit**.

9.5 Connecting probe to MI70 indicator

- ▶ 1. If the probe is installed permanently into a device (for example, an incubator or a chamber), disconnect the probe from the connector.
2. If the MI70 indicator is on, turn it off.
3. Connect the probe to the MI70 indicator using the MI70 connection cable (Vaisala item code: CBL210472).
4. Turn on the MI70 indicator (time and date are requested at first startup). MI70 detects the probe and proceeds to show the measurement screen. The parameters measured by probe will start to show valid measurement results after a few seconds.

9.6 MI70 indicator parts

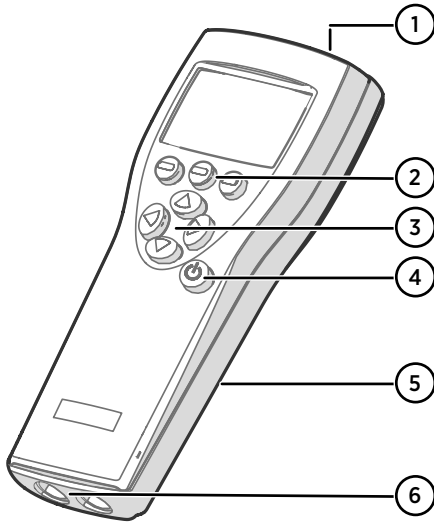







Figure 30 MI70 indicator parts

- 1 Charger socket
- 2 Function key shortcut buttons . The functions change according to what you are doing with the indicator.
- 3 Arrow buttons:
 -  Move up in a menu
 -  Move down in a menu
 -  Enter a sub-menu
 -  Return to previous menu level
- 4 Power On/Off button
- 5 Battery compartment at the back of the indicator
- 6 2 ports with M12 4-pin female connectors for connecting probes and instruments. The ports are labeled I and II.

To open menus, press an arrow button and then press the shortcut buttons. To activate a function shown above the shortcut button, press the shortcut button. To navigate in the menus, press the arrow buttons.

9.7 Holding and saving the display

With the **Hold/Save** function, you can freeze a certain display reading. This reading can be saved in the MI70 memory and it will be available even after MI70 is disconnected from the transmitter.

- ▶ 1. In the basic display, select **Hold/Save**. Alternatively, select **Main menu > Display > Hold/Save display > Hold**.
2. Press **Save**.
3. To view the saved display, go to basic display and select **Record > View recorded data**. Alternatively, select **Main menu > Recording/Viewing > View recorded data**.


A list of saved displays and data recordings appears. The icons on the left of the date and time indicate whether the file is a saved display or a longer recording of data:



Saved display

Data recording

4. Select the saved display based on date and time by pressing the right arrow button.

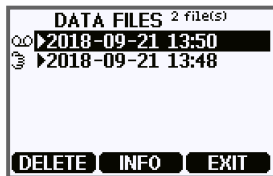

9.8 Recording data

With MI70, you can record transmitter measurement data over a certain period at chosen intervals. These recordings are saved in the MI70 memory and are available even after MI70 is disconnected from the transmitter. To start recording, select the **Record** function key in the basic display, or navigate to the recording menu: **Main menu > Recording/Viewing > Record data**.

9.9 Changing environmental compensation settings with MI70 indicator

You can see the compensation values that are currently used by the probe by selecting them as display quantities from **Main menu > Display > Quantities and Units**. The quantities are as follows:

- **Tcomp**: currently active temperature compensation value.
- **Pcomp**: currently active pressure compensation value.
- **Ocomp**: currently active oxygen concentration compensation value.
- **Hcomp**: currently active relative humidity compensation value.

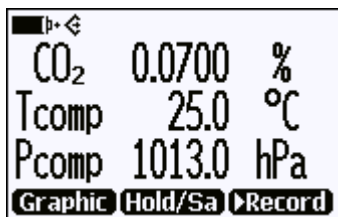


Figure 31 CO₂ reading with Tcomp and Pcomp on MI70 screen

You can change the compensation settings from **Main menu > Settings > Measurement settings**.

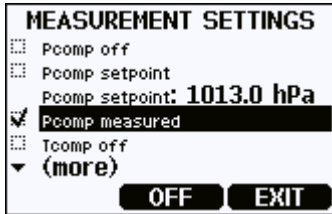


Figure 32 Probe compensation settings on MI70 screen



When you turn a compensation off, the probe still shows a value for the corresponding display quantity (for example, Pcomp shows 1013.2 hPa). This is the default compensation value that is mathematically neutral for the probe's internal compensation model.



GMP251 cannot receive the compensation values from a second probe connected to MI70. Only direct configuration is supported.

10. Vaisala Industrial Protocol

10.1 Vaisala Industrial Protocol overview

The RS-485 line of the probe provides an implementation of the Vaisala Industrial Protocol that can be used for service and configuration use, or for interfacing with the system to which the probe is integrated. The protocol is a plaintext protocol suitable for use both by human operators and automated systems.

10.2 Serial interface settings

Table 10 Default serial interface settings

Property	Description/Value
Bit rate	19200
Parity	None
Data bits	8
Stop bit	1
Flow control	None

10.3 Physical interface

The physical interface is a non-isolated 2-wire interface. The data lines are RS-485 D- and RS-485 D+. Ground is shared with power supply. The connector is a 5-pin male M12.


10.4 Connecting with a computer



- Vaisala USB service cable (item code 242659)
- Computer with:
 - Windows operating system
 - Terminal application (for example PuTTY, available from vaiala.com/software)
 - Free USB port
 - Driver for Vaisala USB service cable installed (available on the cable installation media and at vaiala.com/software)

The steps below describe how to connect to the probe using the PuTTY terminal application for Windows and a USB computer connection cable. Connecting with a computer allows you to configure and troubleshoot your probe using serial line commands.

- ▶ 1. If you have not used the Vaisala USB cable before, install the driver before attempting to use the cable.
2. Connect the USB serial interface cable between your computer and the M12 connector of the probe.
3. Start the PuTTY application.
4. Select **Connection > Serial & USB** and check that the correct COM port is selected in the **Serial or USB line to connect to** field. If you are using the PuTTY terminal application supplied by Vaisala, you can press the **USB Finder** button to open the Vaisala USB Instrument Finder program.
5. Check that the other serial settings are correct for your connection, and change if necessary. **Flow control** should be set to **None** unless you have a reason to change it.

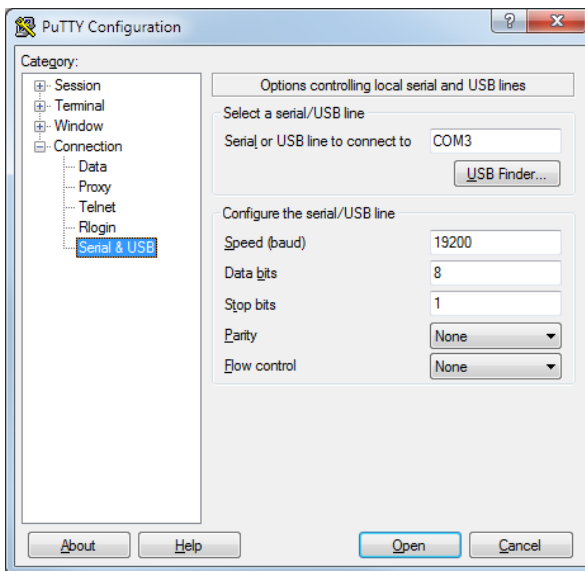

6. Select **Terminal**. Use the following settings:

Local Echo: Select **Force on**. This setting ensures that your typing is shown on the session window.



CAUTION! Make sure that the **Send line ends with line feeds (CR+LF)** option is not selected. Enabling this option can cause issues with the terminal connection.

7. To open the connection window and start using the serial line, select **Open**.




If PuTTY is unable to open the serial port you selected, it shows you an error message instead. If this happens, restart PuTTY and check the settings.

10.4.1 Installing driver for the USB service cable



Only Windows® operating systems are supported by the driver of the USB service cable.

- ▶ 1. Connect the USB service cable to a USB port on your computer. Windows® detects the new device and installs the appropriate driver.
2. Open **Device Manager** from the Windows® Start menu. Use search to find it if necessary (search for "device").
3. Locate the cable in the list of devices:
 - If the device is listed as **Vaisala USB Device** with a COM port number in brackets, the cable is ready for use. **Note the COM port number, you will need it later.**
 - If the device is listed as **Vaisala USB Instrument Cable** without a COM port number listed, you must install the driver manually.

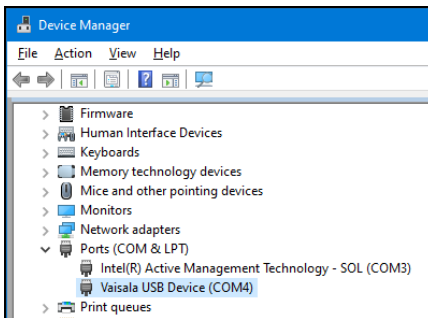

4. To install the driver manually:
 - a. Disconnect the USB service cable from the computer.
 - b. Download the Vaisala USB driver at vaisala.com/software (select the appropriate USB Instrument Driver Setup for your cable).
 - c. Extract the driver zip file and run the USB driver installation program *Vaisala USB Device Driver Setup.exe*. Accept the installation defaults.
 - d. Go back to [step 1](#) and verify that the driver installation works as expected.

10.5 Accessing serial commands from Modbus or analog mode

- ▶ 1. Connect the USB cable to your PC and start the terminal application as instructed in [Connecting with a computer \(page 64\)](#).
2. Start a new terminal session using the default serial settings.
3. Keep the Enter key pressed down and connect the probe to the USB cable.

When the probe is powered on (connected to your PC with the USB cable), you must send five carriage returns (Enter key presses) within 0.7 seconds to force the probe to serial command mode. The probe model information appears in the terminal application when the mode has been successfully changed, and Vaisala Industrial Protocol commands are available for use.
4. To test the connection, enter for example the **?** command. If the mode change failed, close the terminal application, disconnect the probe from the USB cable, and repeat [step 2](#) and [step 3](#).
5. To keep the serial mode in use (forced serial mode access is temporary and switches off at reset), select a serial output option (stop/run/poll) with the **smode** command.



Note that the probe always remains in analog mode when pin #5 is connected to pin #3, and cannot be switched to digital output in this wiring option.

10.6 Enabling Modbus mode from Vaisala Industrial Protocol

If you need to switch from Vaisala Industrial Protocol to Modbus mode, you must configure the following settings:

- Serial line operating mode
- Modbus address
- Serial line settings (bit rate, parity, stop and data bits)

- ▶ 1. Connect the USB cable to your PC and start the terminal application as instructed in [Connecting with a computer \(page 64\)](#).
2. Set the serial mode to Modbus with the **smode** command:

```
smode modbus
```

3. Set the Modbus address to 240 with the **addr** command:

```
addr 240
```

4. Set the serial line settings to 19200/N/8/2 with the **seri** command:

```
seri 19200 N 8 2
```

5. Power off (disconnect) the probe or reset with the **reset** command. The new configuration is available at the next restart.

10.7 Changing from digital output to analog output

- 1. Set up a terminal connection as instructed in [Connecting with a computer \(page 64\)](#).
- 2. Change the mode from digital to analog with the **smode** serial command: **smode analog**.
- 3. Reset the probe (disconnect and reconnect the cable or use the **reset** serial command) to power on in analog output mode.

10.8 Serial commands

The notation **<cr>** refers to the carriage return control character, which you can send in a terminal application by pressing enter on your keyboard. Before entering commands, send a **<cr>** to clear the command buffer.

You can enter the commands in uppercase or lowercase. In the command examples, the keyboard input by the user is in bold type.

[Table 11 \(page 68\)](#) lists the basic serial commands that are available by default. To access advanced serial commands (listed in [Table 12 \(page 69\)](#)), enter the command **pass 1300**.

Table 11 Basic serial commands

Command	Description
Device information and status	
?	Show probe information.
??	Show probe information (will respond in POLL mode).
errs	Show currently active errors.
help	Show list of currently available serial commands.
snum	Show probe serial number.
system	Show probe firmware information.
time	Show probe operation hours and uptime.

Command	Description
vers	Show probe firmware version.
Serial line output and communication	
close	Close connection to probe (POLL mode)
form [modifier string]	Show or set output format.
intv [0 ... 255 s/min/h]	Set continuous output interval for R command.
open [address]	Open connection to probe in POLL mode.
r	Start the continuous outputting.
s	Stop the continuous outputting.
sdelay [1 ... 255]	Show or set serial line transmission delay (1–255, 1 unit = 4 milliseconds)
send	Output a single measurement message.
seri [baud data stop parity]	Show or set the serial interface settings.
smode [mode]	Show or set startup serial mode: RUN , STOP , or POLL .
Environmental compensation	
env	Show or set environmental parameters.
Adjustment information	
adate	Show CO ₂ factory adjustment date.
atext	Show CO ₂ factory adjustment information.
Other commands	
reset	Reset the probe.
pass [1300]	Access advanced serial commands.

Table 12 Advanced serial commands

Command	Description
Serial line output and communication	
addr [0 ... 254]	Show or set probe address.
Analog output	
amode	Show or set analog output mode (analog output limits and error level).
aover	Show or set analog output overrange and clipping behavior.

Command	Description
asel	Show or set analog output parameter and scaling.
Calibration and adjustment	
cco2	Adjust CO ₂ measurement gain and offset.
cdate	Show or set calibration date.
ct	Adjust temperature measurement offset.
ctext	Show or set calibration information.
Environmental compensation	
o2cmode	Show or set oxygen compensation mode.
pcmode	Show or set pressure compensation mode.
rhcmode	Show or set humidity compensation mode.
tcmode	Show or set temperature compensation mode.
Other commands	
frestore	Restore probe to factory settings.

10.9 Device information and status

Table 13 ? command

Syntax	Description
?<cr>	Show listing of device information.
??<cr>	Show listing of device information even if device is in poll mode and connection has not been opened using the open command.
<p>Example:</p> <pre> ? Device : GMP251 Copyright : Copyright (c) Vaisala Oyj 2015. All rights reserved. SW Name : GMP251 SW version : 1.0.0 SNUM : GMP233_5_18 SSNUM : S1234567 CBNUM : c1234567 Calibrated : 20150604 @ Vaisala/R&D Address : 0 Smode : STOP </pre>	

Table 14 errs command

Syntax	Description
errs <cr>	Show active error(s). For a list of possible errors and their remedies, see Error messages (page 111) .
Example (no active errors): <pre> errs NO CRITICAL ERRORS NO ERRORS NO WARNINGS STATUS NORMAL </pre>	

Table 15 help command

Syntax	Description
help <cr>	Show list of currently available serial commands.
Example (showing a list of the basic commands): <pre> help ADATE ADDR ATEXT CLOSE ENV ERRS FORM HELP INTV PASS RR ESET RX SDELAY SEND SENDX SERI SMODE SNUM SYSTEM UNIQID TIME VERS </pre>	

Table 16 snum command

Syntax	Description
snum <cr>	Show serial number of the probe.

Syntax	Description
Example:	
<pre>snum SNUM : M0220028</pre>	

Table 17 system command

Syntax	Description
system<cr>	Show probe firmware information.
Example:	
<pre>system Device Name : GMP251 SW Name : GMP251 SW version : 1.0.0 Operating system : TSF0S1.0</pre>	

Table 18 time command

Syntax	Description
time<cr>	Show how long the probe has been in operation since the last startup or reset. The operation counter is in format hh:mm:ss.
Example:	
<pre>time Time : 01:41:24</pre>	

Table 19 vers command

Syntax	Description
vers<cr>	Show firmware version of the probe.
Example:	
<pre>vers SW version : 1.0.0</pre>	

10.10 Serial line output and communication

Table 20 addr command

Syntax	Description
addr<cr>	Show current device address. Addresses are required for POLL mode.
addr [aaa]<cr>	Set new device address. aaa = address, 0–254 (default = 240)
Example (shows 0 as current address, enter 5 as the new address):	
<pre>addr Address : 0 addr 5 Address : 5</pre>	

Table 21 close command

Syntax	Description
close<cr>	Close the connection that was opened with the open command.
Example:	
<pre>close line closed</pre>	

Table 22 form command

Syntax	Description
form<cr>	Show the currently used measurement format.
form /<cr>	Reset measurement format to default.
form [sss]<cr>	<p>Set a new measurement format.</p> <p>sss = String consisting of modifiers and abbreviations for measured parameters.</p> <p>See Table 23 (page 75) and Table 24 (page 75).</p> <p>Maximum length is 150 characters. Maximum length may be shorter when text strings are used.</p>

Syntax	Description
Example (show currently used measurement format (default format shown here)):	
form 6.0 "C02=" C02 " " U3 #r #n	
Output example (continuous output from RUN mode):	
C02= 452 ppm	
Example (set output format as %CO ₂):	
form 3.1 "C02=" C02% " " U4 #r #n OK	
Output example (continuous output from RUN mode):	
C02= 5.1 %C02 C02= 5.1 %C02 C02= 5.0 %C02 ...	
Example (set output format as CO ₂ ppm with Modulus-65536 checksum):	
form 6.0 "C02=" C02 " " U3 " " CS4 #r #n OK	
Output example (continuous output from RUN mode):	
C02= 3563 ppm 9F C02= 3562 ppm 9E C02= 3559 ppm A4 ...	
Example (set output format as CO ₂ ppm, with start of text (ASCII character 002) and end of text (003) ASCII codes, and without line feed and carriage return at the end):	
form #002 6.0 "C02=" C02 " " U3 #003 OK	
Output example (continuous output from RUN mode, ASCII codes not visible here):	
C02= 866 ppm C02= 866 ppm C02= 867 ppm C02= 867 ppm C02= 867 ppm C02= 868 ppm C02= 868 ppm C02= 869 ppm ...	

Table 23 Output parameters for form command

Output parameter	Abbreviation in form command
Carbon dioxide in ppm	co2
Carbon dioxide in percent	co2%
Currently used temperature compensation value	tcomp
Currently used pressure compensation value	pcomp
Currently used oxygen concentration compensation value	o2comp
Currently used relative humidity compensation value	rhcomp

Table 24 Modifiers for form command

Modifier	Description
x.y	Length modifier (number of digits and decimal places).
#t	Tabulator.
#r	Carriage return.
#n	Line feed.
""	String constant, length 1–15 characters.
#xxx	ASCII code value (decimal) of a special character; for example, #027 for ESC.
addr	Probe address (0–254).
sn	Probe serial number.
ttime	Cumulative operating hours of the probe.
ux	Name of the measurement unit using x number of characters. For example, u3 shows the name of the measurement unit with three characters.
cs4	Modulus-65536 checksum of message sent so far, ASCII encoded hexadecimal notation.
csx	NMEA xor-checksum of message sent so far, ASCII encoded hexadecimal notation.



You can also use the backslash character \ instead of the hash character #.

Table 25 intv command

Syntax	Description
intv<cr>	Show the output interval of the automatically repeating measurement messages (r command and run mode).
intv [iii uu] <cr>	Set the output interval. <ul style="list-style-type: none"> iii = interval, range 0-255 u = unit for interval setting: <ul style="list-style-type: none"> s = seconds min = minutes h = hours If you set the interval to 0, the output messages are output as quickly as they are generated, without additional delay.
Example: <pre>intv 5 s Output interval: 5 S</pre>	

Table 26 open command

Syntax	Description
open [aaa] <cr>	Open a connection to a device at the specified address. Required when device is in poll mode. aaa = address, range 0-254.
Example (target probe in POLL mode, with address 52): <pre>open 52 GMP251: 52 Opened for operator commands</pre>	

Table 27 r command

Syntax	Description
r<cr>	Start the continuous outputting of measurement values as an ASCII text string to the serial line. The probe keeps outputting measurement messages at the interval that has been set with the intv command until stopped with the s command.

Syntax	Description
Example:	
<pre> r C02= 5.1 %C02 C02= 5.1 %C02 C02= 5.1 %C02 C02= 5.0 %C02 C02= 5.0 %C02 ... </pre>	

Table 28 s command

Syntax	Description
s<cr>	Stop the continuous outputting that was started with the r command.
Example:	
<pre> ... C02= 5.1 %C02 C02= 5.0 %C02 C02= 5.0 %C02 s </pre>	

Table 29 sdelay command

Syntax	Description
sdelay<cr>	Show serial line transmission delay (1–255, 1 unit = 4 milliseconds).
sdelay [delay]<cr>	Set a new serial line transmission delay. delay = Serial line delay, range 1–255 (1 unit = 4 milliseconds).
Example (set delay to 0.1 seconds = 25 units):	
<pre> sdelay 25 COM transmit delay : 25 </pre>	

Table 30 send command

Syntax	Description
send<cr>	Output a single measurement message.

Syntax	Description
send [aaa]<cr>	Output a single measurement message from a device in poll mode. aaa = address of the probe, range 0-254
Example: send C02= 5.0 %C02	
Example (target probe in POLL mode, with address 52): send 52 C02= 5.0 %C02	

Table 31 seri command

Syntax	Description
seri<cr>	Show current serial line settings.
seri [b p d s]<cr>	Set new serial line settings. The new settings will be taken into use when the probe is reset or powered up. b = baud rate (9600, 19200, or 38400) p = parity <ul style="list-style-type: none">• n = none• e = even• o = odd d = data bits (7 or 8) s = stop bits (1 or 2) For Modbus, baud rate must be 9600-38400 and parity must be none.
Example (show current settings): seri Com1 Baud rate : 19200 Com1 Parity : N Com1 Data bits : 8 Com1 Stop bits : 1	

Syntax	Description
<p>Example (set serial line to 9600 baud, even, 7 data bits, and 1 stop bit, and reset the probe to take the new settings in use):</p> <pre> seri 9600 e 7 1 OK seri Com1 Baud rate : 9600 Com1 Parity : E Com1 Data bits : 7 Com1 Stop bits : 1 reset GMP251 1.0.0 </pre>	

Table 32 smode command

Syntax	Description
smode<cr>	Show current start-up operating mode of the serial line, and prompt to enter new mode.

Syntax	Description
smode [mode] < cr >	<p>Set serial line start-up operating mode. New mode is taken into use when the device is reset or powered up.</p> <p>Available modes:</p> <p>stop = No automatic output. All commands available. Default mode.</p> <p>run = Automatic output of measurement messages. You can stop the output with the s command, and recontinue with the r command.</p> <p>poll = No automatic output. Will respond to addressed send command and ?? command. You can use other commands after opening a connection using an addressed open command. Use with RS- 485 buses where multiple probes can share the same line.</p> <p>modbus = Serial line communication uses the Modbus protocol (see Modbus (page 31)). Serial line commands (Vaisala Industrial Protocol) are not accessible in the Modbus mode. For instructions on returning to serial mode, see Accessing serial commands from Modbus or analog mode (page 67).</p> <p>analog = Switches the probe from digital output to analog output (active after probe reset). Serial line commands are not accessible in the analog mode: for instructions on returning to serial mode, see Accessing serial commands from Modbus or analog mode (page 67).</p>
<p>Example (set serial mode to "poll"):</p> <div><pre>smode poll Serial mode : POLL</pre></div>	

10.11 Analog output

Table 33 amode command


Syntax	Description
amode [channel]<cr>	<p>Show currently set analog output limits and error level.</p> <p>channel = Analog output channel</p> <ul style="list-style-type: none"> • 1 = voltage output (V) • 2 = current output (mA)
amode [channel lo_value hi_value error_value]<cr>	<p>Set new analog output limits and error output value.</p> <p>channel = Analog output channel</p> <ul style="list-style-type: none"> • 1 = voltage output (V) • 2 = current output (mA) <p>lo_value = Low limit of the channel.</p> <p>hi_value = High limit of the channel.</p> <p>error_value = Error value of the channel.</p> <div>  <p>When setting the error output voltage, note that the probe's maximum voltage output is 10.325 V.</p> </div>
<p>Example (show current configuration):</p> <pre> pass 1300 amode 1 Aout 1 range (V) : 0.00 ... 10.00 (error : 0.00) amode 2 Aout 2 range (mA) : 4.00 ... 20.00 (error : 2.00) </pre>	
<p>Example (set channel 1 to 0–5 V, and error output to 0.0 V; set channel 2 to 0–20 mA, and error output to 23 mA):</p> <pre> amode 1 0 5 0.0 Aout 1 range (V) : 0.00 ... 5.00 (error : 0.00) amode 2 0 20 23 Aout 2 range (mA) : 0.00 ... 20.00 (error : 23.00) </pre>	

Table 34 aover command

Syntax	Description
aover [channel <cr>	Show the behavior of the analog output when the measured value is outside the scaled output range. channel = Analog output channel <ul style="list-style-type: none"> • 1 = voltage output (V) • 2 = current output (mA)
aover [channel clipping error_limit]<cr>	Set the behavior of the analog output when the measured value is outside the scaled output range. channel = Analog output channel <ul style="list-style-type: none"> • 1 = voltage output (V) • 2 = current output (mA) clipping = Output margin (%) at which the output is clipped. error_limit = Measurement value margin (%) at which the output of the channel goes into the error state. The current or voltage output of the error state is defined using the amode command, see Table 33 (page 81) on the previous page.
<p>Example (view currently set analog output overrange behavior on channel 1):</p> <pre> pass 1300 aover 1 Aout 1 clipping :5.00 % Aout 1 error limit :10.00 % </pre>	

Example (for channel 1):

1. View currently set analog output scaling (**asel** command), limits and error level (**amode** command), and overrange behavior (**aover** command):

```

pass 1300
asel 1
Aout 1 quantity : C02(0 ... 200000)
amode 1
Aout 1 range (V) : 0.00 ... 5.00 (error : 0.00)
aover 1
Aout 1 clipping :1.00 %
Aout 1 error limit :5.00 %

```

2. Set analog output overrange clipping to 5 % and error limit to 10 %:

```
aover 1 5 10
Aout 1 clipping : 5.00 %
Aout 1 error limit : 10.00 %
```

The analog output now behaves like this:

- Clipping is now set to 5.00 %, meaning the voltage output is allowed to vary between 0–5.25 V. The analog channel will output the measurement for 0–210 000 ppm, but range 0–5 V remains scaled to 0–200 000 ppm.
- Error limit is 10 %, which means the output will show the error state (0 V) when the measured CO₂ concentration is 10 % outside the scaled output range. With the settings above, this will happen if the measured CO₂ concentration is outside range 0–220 000 ppm.
- The voltage output will never be above 5.25 V because of clipping: the voltage output is clipped when the output reaches 5.25 V, and if the measured CO₂ concentration keeps rising above 220 000 ppm, the output jumps directly to the error state 0 V.

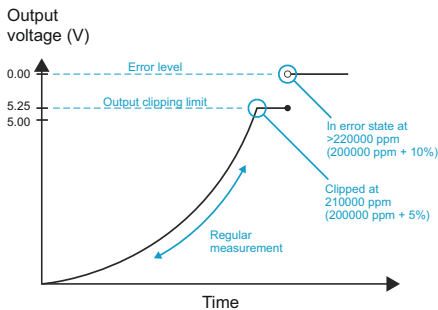


Figure 33 Example of analog output overrange behavior

The example shown above uses output scaled to 0–5 V and 0–200 000 ppm, has error level set to 0 V, clipping set to 5 % overrange, and error limit set to 10 % overrange. CO₂ concentrations (ppm) are indicated for the clipping point and error limit point.

Table 35 asel command

Syntax	Description
asel [channel] <cr>	<p>Show the parameter and scaling of the analog output in ppm.</p> <p>channel = Analog output channel</p> <ul style="list-style-type: none"> • 1 = voltage output (V) • 2 = current output (mA)

Syntax	Description
asel [channel] [parameter lowlimit highlimit] <cr>	Set the parameter and scaling of the analog output. channel = Analog output channel parameter = Parameter that is output on analog channel. The only parameter available is CO ₂ (in ppm). lowlimit = Lower limit of channel scaling in ppm. Minimum value is -1000000 ppm (= -100 %). highlimit = High limit of channel scaling in ppm. Maximum value is 1000000 ppm (= 100 %).
Example (for channel 1, show the currently set analog output parameter and scaling):	
<pre>pass 1300 asel 1 Aout 1 quantity : CO2(0 ... 200000 ppm)</pre>	
Example (for channel 1, set scaling to 0-10 % (= 100 000 ppm)):	
<pre>pass 1300 asel 1 co2 0 100000 Aout 1 quantity : CO2(0 ... 100000 ppm)</pre>	

10.12 Calibration and adjustment commands



CAUTION! Before using the calibration and adjustment commands, read through [Calibration and adjustment \(page 94\)](#). Make sure that the environmental compensation settings of the probe are properly set for your calibration environment; see [Environmental compensation \(page 19\)](#).

Table 36 adate command

Syntax	Description
adate<cr>	Show CO ₂ factory adjustment date.
Example:	
<pre>adate Adjustment date : 20150420</pre>	

Table 37 atext command

Syntax	Description
atext<cr>	Show CO ₂ factory adjustment information.
Example:	
<pre>atext Adjusted at Vaisala/Helsinki</pre>	

Table 38 cdate command

Syntax	Description
cdate<cr>	Show calibration date.
cdate [yyyymmdd]<cr>	Set a new calibration date. yyyymmdd = Year (yyyy), month (mm) and day (dd) of calibration
Example:	
<pre>pass 1300 cdate Calibration date : 20150220</pre>	
Example (set a new calibration date to June 30, 2015):	
<pre>cdate 20150630 Calibration date : 20150630</pre>	

Table 39 ctext command

Syntax	Description
ctext<cr>	Show calibration information text.
ctext [text]<cr>	Set a new calibration information text to be shown after the automatic text "Calibrated at".
Example:	
<pre>pass 1300 ctext Calibrated at 5.0% in lab</pre>	

Syntax	Description
Example (set a new information text):	
<pre>ctest 0% 5% by NN Calibrated at 0% 5% by NN</pre>	

Table 40 cco2 command

Syntax	Description
cco2<cr>	Show current user adjustment status.
cco2 -lo [co2]<cr> cco2 -hi [co2]<cr>	Perform a 1-point (only either low or high concentration) or 2-point (both low and high concentrations) calibration and adjustment. -lo = Adjustment at low concentration (under 2 %CO2) -hi = Adjustment at high concentration (over 2 %CO2) co2 = CO2 concentration reference in ppm
cco2 -save<cr>	Save the currently entered adjustments. Successfully saving the adjustment clears the calibration date (cdate command) and calibration text (ctest command) that have been stored in the probe. Use those commands to enter a new calibration date and text.
cco2 -cancel<cr>	Cancel currently entered adjustments.
cco2 -reset<cr>	Clear user adjustments.
Example (show current user adjustment status; no adjustment done):	
<pre>pass 1300 cco2 1.Ref. point low 0 1.Meas. point low 0 2.Ref. point high 200000 2.Meas. point high 200000 Gain : 1.0000 Offset : 0.0000</pre>	

Syntax	Description
<p>Example (perform a 1-point calibration):</p> <ol style="list-style-type: none"> 1. Let the probe stabilize in the desired CO₂ concentration (here: 5 %CO₂ (=50000 ppm)). 2. Enter the calibration commands: <pre>pass 1300 cco2 -hi 50000 OK cco2 -save OK</pre> <ol style="list-style-type: none"> 3. Enter a new calibration date and information text: <pre>cdate 20150630 Calibration date : 20150630 ctext 5% in lab Calibrated at 5% in lab</pre>	
<p>Example (perform 2-point calibration):</p> <ol style="list-style-type: none"> 1. Let the probe stabilize in the desired low CO₂ concentration (here: 0 %CO₂). 2. Enter the calibration commands: <pre>pass 1300 cco2 -lo 0 OK cco2 -save OK</pre> <ol style="list-style-type: none"> 3. Let the probe stabilize in the desired high CO₂ concentration (here: 5 %CO₂ (=50000 ppm)). 4. Enter the calibration commands: <pre>pass 1300 cco2 -hi 50000 OK cco2 -save OK</pre> <ol style="list-style-type: none"> 5. Enter a new calibration date and information text: <pre>pass 1300 cdate 20150430 Calibration date : 20150430 ctext 0% 5% by NN Calibrated at 0% 5% by NN</pre>	

10.13 Environmental compensation commands




To apply an accurate relative humidity compensation, the temperature and pressure compensation configurations must also match your measurement environment. See the **rhcmode**, **tcmode** and **pcmode** commands for instructions on enabling compensation configuration, and **env** command for instructions on setting a compensation value.



For more information on environmental compensation and the default (neutral) compensation values used for disabled compensations, see [Environmental compensation \(page 19\)](#).

Table 41 env command

Syntax	Description
env<cr>	Show current compensation values. Before using this command, you must enable environmental compensation using the following commands: <ul style="list-style-type: none">• o2cmode [on]• pcmode [on]• rhcmode [on]• tcmode [on measured]
env [temp pres oxy hum] [value]<cr>	Set new permanent compensation values and store them in EEPROM . EEPROM: <ul style="list-style-type: none">• Non-volatile memory, values retained during power off.• Number of writes is limited to 30000 cycles by memory implementation.• Must only be used for writing permanent values, to avoid wearing out the EEPROM. temp = Compensation temperature. Range -40 ... +100 °C (-40 ... +212 °F). pres = Compensation pressure. Range 500–1100 hPa. oxy = Oxygen content of background gas. Range 0–100 %. hum = Relative humidity of background gas. Range 0–100 %.

Syntax	Description
env [xtemp xpres xoxy xhum] [value] <cr>	<p>Set new compensation values and store them in RAM.</p> <p>RAM:</p> <ul style="list-style-type: none"> • Volatile memory that loses the values when probe is reset, and where values are loaded from non-volatile (EEPROM) memory at power-on. • Must be used for continuously updated values. <p>xtemp = Compensation temperature stored in RAM. Range -40 ... +100 °C (-40 ... +212 °F).</p> <p>xpres = Compensation pressure stored in RAM. Range 500-1100 hPa.</p> <p>xoxy = Oxygen content of background gas stored in RAM. Range 0-100 %.</p> <p>xhum = Relative humidity of background gas stored in RAM. Range 0-100 %.</p> <div>  <p>Note: If temperature compensation is configured to use an internally measured value (tcmode is set to measured), it will continuously update the value in RAM, overriding any temperature value that is written to RAM with the ENV command.</p> </div>
<p>Example (Show current compensation values; all compensations are enabled. Note that temperature compensation is in "measured" mode, so the value in use is constantly changing):</p> <pre> env In eeprom: Temperature (C) : 25.00 Pressure (hPa) : 1013.00 Oxygen (%O2) : 21.00 Humidity (%RH) : 50.00 In use: Temperature (C) : 36.40 Pressure (hPa) : 1013.00 Oxygen (%O2) : 19.70 Humidity (%RH) : 93.00 </pre>	

Syntax	Description
Example (set temperature compensation to setpoint mode, and change temperature setpoint value to 37.2 in RAM):	
<pre>pass 1300 tcmode on T COMP MODE : ON env xtemp 37.2 In eeprom: Temperature (C) : 25.00 Pressure (hPa) : 1013.00 Oxygen (%O2) : 21.00 Humidity (%RH) : 50.000 In use: Temperature (C) : 37.2 Pressure (hPa) : 1013.00 Oxygen (%O2) : 19.70 Humidity (%RH) : 93.00</pre>	


Table 42 o2cmode command

Syntax	Description
o2cmode<cr>	Check current oxygen compensation mode. Possible modes: <ul style="list-style-type: none">• on = Compensation enabled using setpoint value.• off = Compensation disabled, default (neutral) value used: see Environmental compensation (page 19)
o2cmode [on off]<cr>	Change oxygen compensation mode (on or off).
Example (check oxygen compensation mode; oxygen compensation is disabled, a neutral value is used):	
<pre>pass 1300 o2cmode O2 COMP MODE : OFF</pre>	
Example (enable oxygen compensation):	
<pre>pass 1300 o2cmode on O2 COMP MODE : ON</pre>	

Table 43 pcmode command

Syntax	Description
pcmode<cr>	Check current pressure compensation mode. Possible modes: <ul style="list-style-type: none"> • on = Compensation enabled using setpoint value. • off = Compensation disabled, default (neutral) value used: see Environmental compensation (page 19).
pcmode [on off]<cr>	Change pressure compensation mode (on or off).
Example (check pressure compensation mode; pressure compensation is enabled using a setpoint value):	
<pre>pass 1300 pcmode P COMP MODE : ON</pre>	

Table 44 rhcmode command

Syntax	Description
rhcmode<cr>	Check current relative humidity compensation mode. Possible modes: <ul style="list-style-type: none"> • on = Compensation enabled using setpoint value. • off = Compensation disabled, default (neutral) value used: see Environmental compensation (page 19). <div>  To apply an accurate relative humidity compensation, make sure that the temperature compensation and pressure compensation configurations also match the measurement environment. </div>
rhcmode [on off]<cr>	Change relative humidity compensation mode (on or off).

Syntax	Description
Example (check relative humidity compensation mode; relative humidity compensation is disabled, a neutral value is used):	
<pre>pass 1300 rhcmode RH COMP MODE : OFF</pre>	
Example (enable temperature, pressure and relative humidity compensation using setpoint values):	
<pre>pass 1300 tcmode on T COMP MODE : ON pcmode on P COMP MODE : ON rhcmode on RH COMP MODE : ON</pre>	

Table 45 tcmode command

Syntax	Description
tcmode<cr>	Check current temperature compensation mode. Possible modes: <ul style="list-style-type: none"> • on = Compensation enabled using setpoint value. • off = Compensation disabled, default (neutral) value used: see Environmental compensation (page 19). • measured = Compensation enabled using internal measurement.
tcmode [on off measured]<cr>	Change temperature compensation mode (on , off or measured).
Example (check temperature compensation mode; temperature compensation is enabled using a setpoint value):	
<pre>pass 1300 tcmode T COMP MODE : ON</pre>	
Example (change temperature compensation to use internal measurement):	
<pre>pass 1300 tcmode measured T COMP MODE : MEASURED</pre>	

10.14 Other commands

Table 46 frestore command


Syntax	Description
frestore<cr>	<p>Restore the probe to its factory configuration. All user settings and user calibration parameters will be lost.</p> <div>  <p>After using the frestore command, reset the probe using the reset command.</p> </div>
<p>Example (restore the factory settings and reset the probe):</p> <pre>pass 1300 frestore Parameters restored to factory defaults reset GMP251 1.0.0</pre>	

Table 47 pass command

Syntax	Description
pass [code]<cr>	<p>Access advanced serial commands. Advanced commands can be used until the next reset.</p> <p>code = Code for enabling advanced commands (1300).</p>
<p>Example:</p> <pre>pass 1300</pre>	

Table 48 reset command

Syntax	Description
reset<cr>	Reset the probe. The probe will restart as if it had just been powered on.
<p>Example:</p> <pre>reset GMP251 1.0.0</pre>	

11. Calibration and adjustment

Calibrate and adjust the CO₂ measurement of the GMP251 probe as needed. Before starting, read through this section completely so that you are aware of your options, and the main factors that affect the result.



Performing an accurate calibration and adjustment takes some time and preparation. Instead of doing it yourself, you can also have a Vaisala Service Center calibrate and adjust your probe.



Calibration means comparing the measurement output of the device to a known reference, such as a known environment in a calibration chamber or the output of a reference instrument. Correcting the reading of the device so that it measures accurately is referred to as **adjustment**.

11.1 Calibration setup

Using calibrated probe as reference

You can perform a 1- or 2-point calibration using the Indigo80 indicator and a calibrated reference instrument to compare against, for example a calibrated GMP251 probe connected to Indigo80.

With the probe being calibrated and the reference instrument in the same space, allow the measurements to stabilize before comparing the readings. Try to provide as stable an environment as you can during this time. Avoid working around the probe and reference instrument during this time.

Using calibration gas as reference

There are two easy ways to use a calibration gas as a reference:

- You can supply the gas to the probe using the calibration adapter accessory (Vaisala item code DRW244827SP). Gas flow should be in the range 0.5 – 1.1 l/min, recommendation is 1 l/min. Allow the measurement to stabilize for 3 minutes before starting the calibration.
- You can fill the entire incubator with the calibration gas. You can use calibration gas as a reference by putting the probe in a suitable chamber (for example, an incubator) and filling that chamber with the calibration gas.

To perform a 2-point adjustment, you need 2 calibration gases: one gas that is below 2 %CO₂ (low-end reference) and one that is above 2 %CO₂ (high-end reference).

When supplying the gas from a gas bottle, make sure the gas bottle has stabilized to room temperature before starting.

11.2 Effect of environmental compensations

The probe has various environmental compensations that improve its CO₂ measurement accuracy (see [Environmental compensation \(page 19\)](#)). As the calibration and adjustment environment may differ from the actual measurement environment, you must make sure that the compensation settings are properly set. Here are some key points to remember:

- Pressure and temperature compensations have a significant effect on accuracy. If you are using setpoint values instead of the values from the built-in temperature sensor or an integrated system, make sure to correct the setpoints so that they correspond to your calibration situation. Consider switching temperature compensation to use the internal sensor and/or integrated system when calibrating, and then switching back when calibration and adjustment is done.
- The effect of background gas compensations for humidity and oxygen may be significant when using calibration gases, since these gases are often dry and oxygen-free. For example, pure nitrogen gas is typically used as a convenient 0 ppm CO₂ reference. As it does not contain any oxygen or humidity, the compensations for them must be set to zero.
- Remember to restore the normal compensation settings after completing calibration and adjustment. If you are integrating the calibration functionality of the probe as part of a control software, also implement proper handling of the environmental compensations.

11.3 Limits of adjustment

The probe limits the amount of adjustment that is allowed to the CO₂ measurement. The maximum correction that you can apply is 1000 ppm + 25 % of the probe's uncorrected reading. Previous user adjustments do not affect this limit (the correction is not cumulative). This feature limits the possible error introduced by incorrect adjustment.

For example, if you are adjusting using a 5 %CO₂ calibration gas (50000 ppm), the maximum correction you can apply is approximately 13500 ppm. Attempting to apply a greater adjustment will fail. Notification of failure from the probe depends on the interface you are using for adjustment.

11.4 Adjustment types

You can adjust the CO₂ measurement of the probe in 1 or 2 points.

- 1-point adjustment is recommended if you are interested in maintaining a fixed CO₂ level. For best result, use a calibration gas with a CO₂ concentration that is close to the intended level.
- 2-point adjustment is recommended if you typically measure a variable CO₂ level.

Available adjustment functions depend on the interface you use to operate the probe. If you want to compare the reading of the probe to a reference instrument and adjust it accordingly, use the Indigo80 handheld indicator and a reference probe. If you want to integrate the functionality into a control system, the Modbus interface and the Vaisala industrial protocol are recommended.

Vaisala Insight PC software

You can carry out 1-point and 2-point adjustments and change the environmental compensation settings with the Insight software.

Vaisala Indigo transmitters and Indigo80 handheld indicator

You can use an Indigo transmitter's configuration interface in the Insight PC software or connect the probe to the Indigo80 handheld indicator to carry out 1-point and 2-point adjustments and change the environmental compensation settings.

MI70 handheld indicator

The MI70 handheld indicator supports 1-point adjustment, either using a calibration gas or using a reference instrument that is connected to the MI70.

Vaisala Industrial Protocol

Vaisala industrial protocol supports 1-point and 2-point adjustment with the **cco2** command. Configuration of the environmental compensation settings can be done using serial line commands.

More information

- [Calibration and adjustment commands \(page 84\)](#)
- [Overview of calibration with Indigo200 using Insight \(page 98\)](#)
- [Performing manual adjustment with Indigo80 using predetermined values \(page 103\)](#)
- [Calibration and adjustment with MI70 indicator \(page 104\)](#)

11.5 DRW244827SP calibration adapter

The optional calibration adapter accessory can be used to feed a reference gas to the probe through a gas port when calibrating. Gas flow should be in the range 0.5–1.1 l/min, recommendation is 1 l/min. Allow the measurement to stabilize for 3 minutes before starting the calibration.

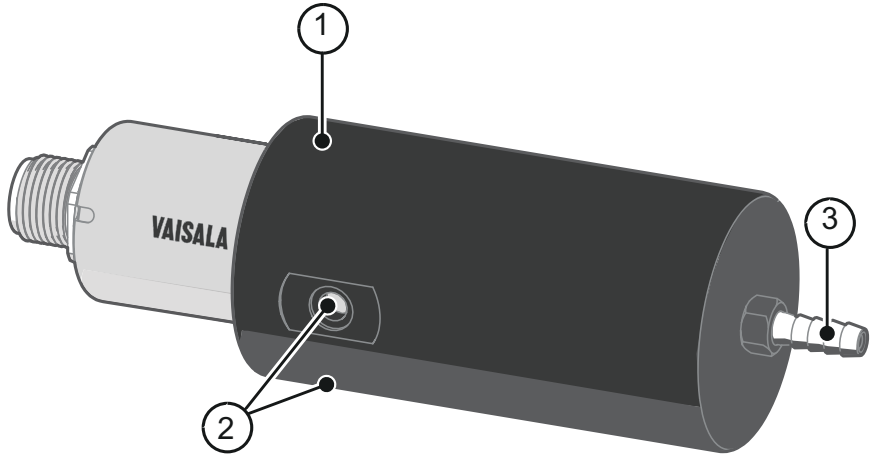


Figure 34 DRW244827SP calibration adapter with probe inserted

- 1 O-ring inside the adapter
- 2 Gas outlet on each side of the adapter
- 3 Gas port (port outer diameter 4.6 mm, port hole inner diameter 2 mm, suitable for tubing with 4 mm inner diameter)

11.6 Overview of calibration with Indigo200 using Insight

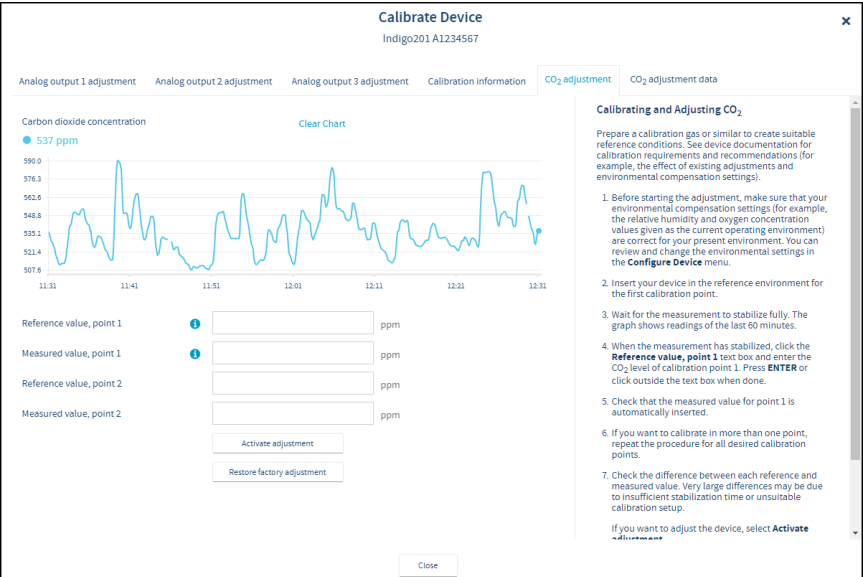


Figure 35 Insight calibration menu main view (CO₂ probe)

The calibration menu in the Insight PC software contains options for calibrating and adjusting the measurement of the probe you have connected to the transmitter, and testing and adjusting the analog output channels of the transmitter. You can also view the current adjustment and restore the probe's factory adjustments.




The range of available options (for example, environmental compensations and the number of adjustment points) varies depending on the features of the connected probe.
The menu examples presented here are examples of calibration options that are available when a CO₂ probe is connected to the transmitter.



CAUTION! Before adjusting a probe's measurement, make sure that you have familiarized yourself with the probe-specific calibration requirements such as possible adjustment limits and environmental compensation interdependencies. See the user documentation of the probe for probe-specific information.

Starting calibration

To open the calibration menu, select  to access Insight main menu. In the main menu, select **Calibrate**.

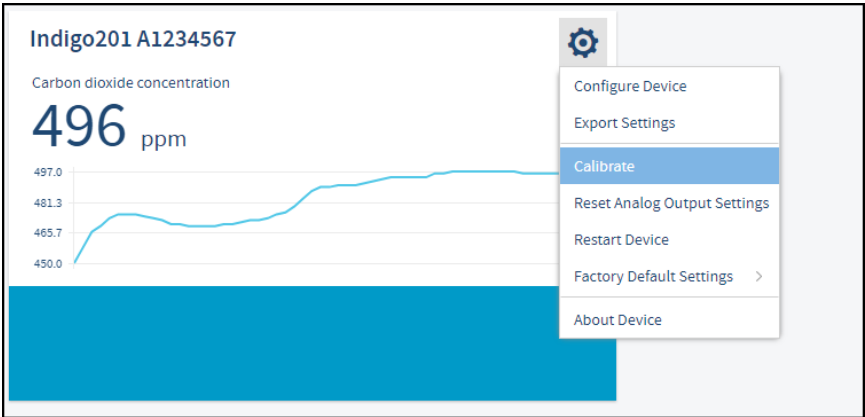


Figure 36 Insight **Calibrate** menu selection

Restoring factory adjustment

Before adjusting the probe measurement, you can remove any possible earlier adjustments with the **Restore factory adjustment** selection, if needed.

The screenshot shows the calibration screen with four input fields for reference and measured values at two points. Each field has an information icon (i) to its left. Below the input fields are two buttons: 'Activate adjustment' and 'Restore factory adjustment'. The 'Restore factory adjustment' button is highlighted with a red rectangle.

Figure 37 Insight **Restore factory adjustment** selection

11.6.1 Environmental compensation

When making adjustments, you can enter information about the environment in which you perform the adjustment. This allows compensating for environmental factors that have an effect on the measurement (for example, temperature, pressure, or background gases).

The selection of environmental compensations available in the Indigo transmitter configuration menus depends on the features and configuration of the connected probe. If the probe you have connected to the transmitter supports environmental compensations, the configuration options are shown in the **Configure device** menu accessible from the Insight PC software main view. For information and instructions on the environmental compensation features of the probe, see the user documentation of the probe you are using.

Configure Device In Advanced Mode

Indigo201 A1234567

Analog output 1

Analog output 2

Analog output 3

Relay output A

Relay output B

Compensation power-up defaults

Compensation setpoints

Measurement

Pressure compensation on/off ☒

Temperature compensation mode Measured ▾

Humidity compensation on/off ☒

Oxygen compensation on/off ☒

Filtering factor

Figure 38 Insight environmental compensation configuration example with CO₂ probe

11.6.2 CO₂ probe adjustment example

Prepare a calibration gas or similar to create suitable reference conditions. See device documentation for calibration requirements and recommendations (for example, the effect of existing adjustments and environmental compensation settings).

To adjust the CO₂ measurement of the probe that is connected to the transmitter:

- ▶ 1. Before starting the adjustment, make sure that your environmental compensation settings (for example, the relative humidity and oxygen concentration values given as the current operating environment) are correct for your present environment. You can review and change the environmental settings in the **Configure device** menu.
2. To open the Insight PC software calibration menu, select to access Insight main menu. In the main menu, select **Calibrate**.
3. Select **CO₂ adjustment**.
4. Insert the probe in the reference environment for the first calibration point.
5. Wait for the measurement to stabilize fully. The graph shows readings of the last 60 minutes.

6. When the measurement has stabilized, click the **Reference value, point 1** text box and enter the CO₂ level of calibration point 1. Press **ENTER** or click outside the text box when done.
7. Check that the measured value for point 1 is automatically inserted.
8. If you want to calibrate in more than one point, repeat the procedure for all desired calibration points.
9. Check the difference between each reference and measured value. Very large differences may be due to insufficient stabilization time or unsuitable calibration setup.

If you want to adjust the device, select **Activate adjustment**.

To exit without taking the adjustment in use, select **Close**.

10. After calibrating your device, update the information in the **Calibration information** tab.

11.7 Calibration and adjustment with Indigo80 handheld indicator

The Indigo80 handheld indicator provides two ways to calibrate and adjust measurements for probes: **Guided calibration** and **Manual adjustment**.

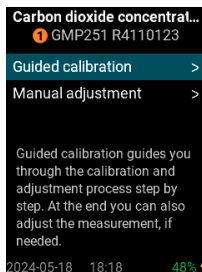


Figure 39 Calibration methods available with Indigo80

Use **Guided calibration** if you have a reference gas environment or a calibrated reference instrument (another GMP251 or GMP252 probe or a third party instrument), and you want to:

- Perform calibration check at 1 or 2 CO₂ levels (check only, probe not adjusted), **or**
- Adjust your probe using 1 or 2 CO₂ levels (probe adjusted to minimize the observed error).

Use the **Manual adjustment** mode only if you already have calibration data (for example, a calibration certificate) for 1 or 2 calibration points and want to adjust the probe to correct the observed errors. No reference instrument or calibration gases are needed. For all other use cases, **Guided calibration** is recommended.

11.7.1 Performing guided calibration with Indigo80



- Probe connection cable
- Calibration reference (for example, a reference gas with a known concentration, or a recently calibrated Indigo80-compatible probe)

In guided calibration, the indicator guides you through the calibration process step by step. At the end, you can also adjust the measurement, if needed.

You can choose to perform either 1-point or 2-point calibration, and choose to use the reference gas or instrument of your choice, for example, another probe connected to your indicator.

Before starting the adjustment:

- Familiarize yourself with information on adjustment points and requirements for your probe.
- Make sure that the probe's environmental compensation settings (for example, the pressure value given as the current operating pressure) are correct for your present environment.
- **Always restore factory adjustment before entering a new adjustment.** This prevents any possible earlier adjustments having an effect on the new adjustment you make.

1. Connect the probe(s) to the indicator.
2. Open the indicator main menu by pressing
3. Select **Devices**, then select the probe to be calibrated.
4. Select **Calibrate / Adjust**.
5. Read the notification and select **OK** to close it.
6. Select **CO₂ adjustment**.
7. Select **Guided calibration**.
8. In **Calibration points**, select the number of calibration points: **1 point** or **2 points**.
9. In **Reference**, select the reference instrument.
Selection **Other reference** allows you to use any reference environment(s) of your choice (for example, a reference gas).
10. Select **Start calibration**.

The indicator guides you through the calibration process from this point onwards. Press to proceed from one view to another.

More information

- [Configuring environment settings with Indigo80 handheld indicator \(page 48\)](#)
- [Effect of environmental compensations \(page 95\)](#)
- [Restoring factory adjustment with Indigo80 handheld indicator \(page 104\)](#)

11.7.2 Performing manual adjustment with Indigo80 using predetermined values



- Probe connection cable

This section tells you how to use the Indigo80 indicator to adjust the probe's measurement in 1 or 2 points. Use the **Manual adjustment** mode to enter predetermined calibration data (for example, from a calibration certificate) to adjust your probe to correct the observed errors.

- ▶ 1. Connect the probe to the indicator.
2. Open the indicator main menu by pressing \ominus .
3. Select **Devices**, then select the probe to be calibrated.
4. Select **Calibrate / Adjust**.
5. Read the notification and select **OK** to close it.
6. Select **CO₂ adjustment**.
7. Select **Manual adjustment**. Read the notification and select **OK** to close it.
8. Select **Enter manually**.
 - a. Select **Reference value, point 1**.
 - b. Select **ppm** or **%** as the unit.
 - c. Use the arrow buttons to set the reference value, then select **Enter** or **Set**.
 - d. In the **CO₂ adjustment** view, the probe prefills the **Measured value, point 1** field. Select the field and change the value to match your actual predetermined value (the value that this probe reported in the given reference environment).

Scroll down and repeat for **Reference value, point 2** and **Measured value, point 2** (if performing a 2-point adjustment).
9. To store the adjustments to the probe, scroll down and select **Activate adjustment > OK**. Confirm by selecting **OK**.

To exit without taking the adjustment in use, press \ominus in the **CO₂ adjustment** view to end calibration.
10. In the **Calibrate / Adjust** menu of your probe, select **Calibration information** and update any relevant data.
11. Press \ominus to end calibration.

More information

- [Configuring environment settings with Indigo80 handheld indicator \(page 48\)](#)
- [Effect of environmental compensations \(page 95\)](#)
- [Restoring factory adjustment with Indigo80 handheld indicator \(page 104\)](#)

11.7.3 Restoring factory adjustment with Indigo80 handheld indicator



- Probe connection cable

If needed, you can restore the factory adjustment for your probe's output parameters using the Indigo80 handheld indicator.



Always restore factory adjustment before entering a new adjustment. This prevents any possible earlier adjustments having an effect on the new adjustment you make.

1. Connect the probe(s) to the indicator.
2. Open the indicator main menu by pressing
3. Select **Devices**. If you have more than one device connected to the indicator, make a further selection between the devices.
4. Select **Calibrate / Adjust**.
5. Select the type of adjustment to restore.
For example, select **CO₂ adjustment**.
6. Select **Manual adjustment**.
7. Select **Enter manually**.
8. Scroll down and select **Restore factory adjustment > OK**.
9. Confirm by selecting **OK**.
10. Press to end calibration.

11.8 Calibration and adjustment with MI70 indicator




Before using the MI70 indicator for calibration and adjustment, read the instructions in [Calibration and adjustment \(page 94\)](#). Make sure that the environmental compensation settings of the probe are properly set for your calibration environment; see [Changing environmental compensation settings with MI70 indicator \(page 62\)](#).



When two probes are connected to the MI70 indicator, MI70 uses Roman numerals “I” and “II” to indicate which port the parameter or function in question is connected to.

11.8.1 1-point adjustment with an MI70-compatible reference probe

- ▶ 1. Connect the GMP251 probe to Port I of the MI70 indicator.
2. Connect the calibrated reference probe to Port II. Make sure the reference probe is in the same environment as the GMP251's sensor.
3. If you are using the calibration adapter to feed a calibration gas to the GMP251, you must feed the same gas to the reference probe also. Refer to the documentation of your reference probe on how to do this, and what accessories you need.
4. Turn on the MI70 indicator.
5. Start the adjustment sequence from **Main menu > Functions > Adjustments**.

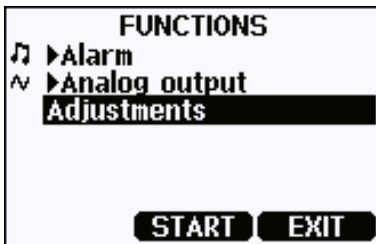


6. MI70 notifies you that automatic power off is disabled during adjustment mode, press **OK** to acknowledge.
7. To proceed with the adjustment, select the **CO2(I)** parameter in the **Select Quantity** screen.


In the **Select Quantity** screen you can also view the currently used compensation values, and the **Last adjustment date** information. You can update the date and text using the **CDATE** and **CTEXT** commands on the serial line.

8. You may be prompted to check the environmental settings of the reference probe before proceeding. Press **Yes** to check the settings and **Exit** when you have checked and corrected the settings.

9. The adjustment mode is now active, and you can see the measured CO₂ readings and their difference on the screen. Allow the measurement to stabilize. To proceed with the adjustment, press **Adjust**.

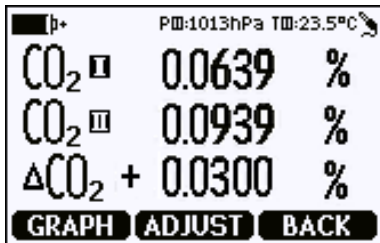

10. Select **To same as CO₂(II)**.

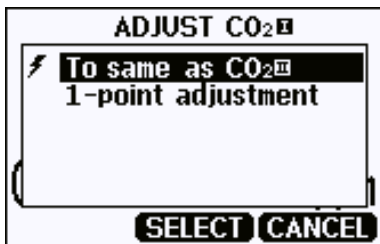

11. You will be prompted to confirm you want to adjust: select **Yes**.
12. If the adjustment is successful, MI70 will show the text **Adjustment Done**, after which you will return to the adjustment mode. At this point you can press **Back** and **Exit** to leave the adjustment mode. The adjustment is now completed.

If the adjustment cannot be applied, MI70 will show the text **Cannot adjust**, possibly followed by a text stating the reason. A possible reason for an adjustment failure is attempting to apply a very large correction to the reading.

11.8.2 1-point adjustment with a reference gas

- ▶ 1. Connect the GMP251 to Port I of the MI70 indicator.
2. Feed a calibration gas to the GMP251 using the calibration adapter accessory (Vaisala item code: DRW244827SP). If you are using ambient air as the calibration gas, you must have a reference meter in the same environment to verify the CO₂ concentration.
3. Turn on the MI70 indicator.


4. Start the adjustment sequence from **Main menu > Functions > Adjustments**.




5. MI70 notifies you that automatic power off is disabled during adjustment mode, press **OK** to acknowledge.
6. Select the CO₂ parameter when prompted.
7. You may be prompted to check the environmental settings of the reference probe before proceeding. Press **Yes** to check the settings and **Exit** when you have checked and corrected the settings.
8. The adjustment mode is now active, and you can see the measured CO₂ reading on the screen. To proceed with the adjustment, press **Adjust**.



9. Select **1-point adjustment**.

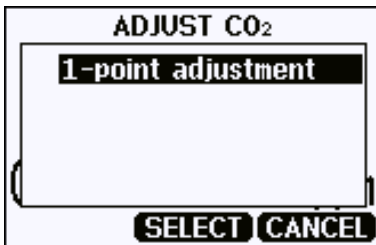

10. You will be prompted if you really want to adjust. Select **Yes**.
11. You are now in the 1-point adjustment screen. Allow the measurement to stabilize and press **Ready**.
12. Enter the CO₂ concentration of the reference gas and press **OK**.

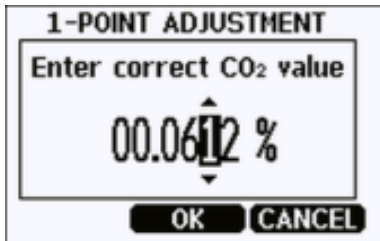

13. You will be prompted if you really want to adjust. Select **Yes**.
14. If the adjustment is successful, MI70 will show the text **Adjustment Done**, after which you will return to the adjustment mode. At this point you can press **Back** and **Exit** to leave the adjustment mode. The adjustment is now completed.

If the adjustment cannot be applied, MI70 will show the text **Cannot adjust**, possibly followed by a text stating the reason. A possible reason for an adjustment failure is attempting to apply a very large correction to the reading.

12. Maintenance

12.1 Cleaning GMP251 and GMP80P

You can clean the probe body by wiping it with a moist cloth. Standard cleaning agents can be used.

When cleaning, follow these precautions:

- Do not immerse the probe in liquid to clean it.
- GMP251 only:
 - Be careful not to block the filter when cleaning the probe. The optional sintered PTFE filter is especially sensitive to blockage.
 - When changing the filter, you can use clean instrument air to gently blow any loose dirt and filter material from the sensor. Do not attempt to clean the optical surfaces in any other manner.

12.1.1 Chemical tolerance



Detach the GMP251 probe from GMP80P before cleaning GMP251 with hydrogen peroxide, acetone, or acetic acid.

The following chemicals can be used to clean the GMP251 probe:

- H_2O_2 (2000 ppm), non-condensing
- Alcohol-based cleaning agents such as ethanol and IPA (70 % isopropyl alcohol, 30 % water)
- Acetone
- Acetic acid

The following chemicals can be used to clean the GMP80P body:

- Alcohol-based cleaning agents such as ethanol and IPA (70 % isopropyl alcohol, 30 % water).



Avoid exposing the probe to chemicals for unnecessarily long periods of time. Do not immerse the probe in a chemical or spray anything directly on the probe. Wipe chemicals off the probe after cleaning.

12.2 Detaching GMP251 from GMP80P


- ▶ 1. Stop the measurement by switching off Indigo80 or by disconnecting GMP80P from Indigo80.
2. Loosen the locking ring of the cable connector on the bottom of GMP251 and disconnect the cable from the probe.

3. Grip the GMP80P body with both hands and push the GMP251 probe out using your thumbs.



Avoid pulling the GMP251 probe out by the cable connector. Pulling and/or twisting the connector may break it.

12.3 Changing the GMP251 filter

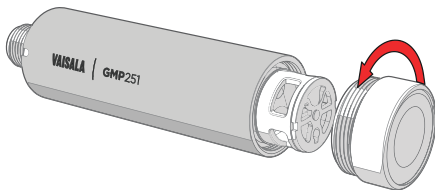

Change the filter to a new one if it shows visible signs of contamination or dirt. When changing the filter, use clean gloves to avoid blocking the pores of the new filter.

Figure 40 Opening the filter



CAUTION! When changing the filter, you can use clean instrument air to gently blow any loose dirt and filter material from the sensor. Do not attempt to clean the optical surfaces in any other manner.

13. Troubleshooting

13.1 Problem situations with GMP251

Problem	Possible cause	Remedy
Analog output reading is unchanging and appears incorrect.	Analog output is in error state.	Remove the cause of the error state and the analog output will recover its normal function.
Probe outputs stars "*****" on serial line instead of measurement data.	Incorrect supply voltage.	Check the power supply. Check the active errors using the ERRS command on the serial line.
	Unsuitable operating environment.	Verify that the operating environment is within specified operating limits.
Unable to access probe on the RS-485 line.	Incorrect wiring.	Check that the RS-485 connection is wired correctly.
	Probe in POLL mode with unknown address.	Power cycle or reset the probe and try again.
CO ₂ measurement not working.	Condensation on the sensor.	Remove the filter and check if condensation has formed on the sensor. If yes, dry out the condensation with instrument air and insert a new dry filter. Keep the probe powered and operating to prevent re-occurrence.

More information

- [Wiring \(page 30\)](#)
- [Analog output error state \(page 113\)](#)

13.2 Error messages

The error messages are categorized according to the severity of the status:

- **Critical errors** are fatal to the operation of the unit. It may not be able to respond to communication at all, and will not measure correctly.
- **Errors** prevent CO₂ measurement and cause the analog outputs to be set to the error state. Depending on the problem, errors may resolve themselves. For example, sensor heating will eventually dry out condensation on the optical surfaces.
- **Warnings** do not prevent normal operation but may indicate possible problems.

- **Status** indicates a known state of the unit.

In case of constant error, contact Vaisala.

Error message and code	Description	Recommended action
Critical errors		
Program memory crc critical error [1]	Program memory is corrupted.	Fatal error, contact Vaisala technical support.
Parameter memory crc critical error [2]	Parameter memory is corrupted.	Fatal error, contact Vaisala technical support.
Errors		
Low supply voltage error [5]		Check and correct the supply voltage.
Internal 30V error [6]	Low internal 30 V voltage.	
Low RX signal error [7]	Low input signal. Can be caused by dirt or condensation on the optical surfaces.	Wait to see if condensation is removed by heat.
Internal 8V error [8]	Low internal 8 V voltage.	
RX signal cut error [9]	Signal distortion (EMC interference)	
Out of measurement range error [13]	CO ₂ concentration is too high to measure.	Wait for CO ₂ concentration to fall into the measurable range.
Sensor heater error [14]	Sensor heater resistance is out of range.	
IR temperature error [15]	Infrared source temperature error.	
FPI slope error [16]	Signal receiver error.	Contact Vaisala technical support.
Internal 2.5V error [17]	Internal 2.5 V voltage out of range.	
Internal 1.7V error [18]	Internal 1.7 V voltage out of range.	
Low IR current error [19]	Infrared source failure.	Contact Vaisala technical support.
Warnings		
Signal too low warning [21]	Low input signal. Can be caused by dirt or condensation on the optical surfaces.	Continue normally.

Error message and code	Description	Recommended action
Cut warning [23]	EMC interference error limit approaching.	Check for EMC interference sources.
Unexpected restart detected [24]	Transmitter is reset by watchdog process.	Continue normally.
Status messages		
CO ₂ adjustment mode active [27]		Complete the CO ₂ adjustment.
Calibration about to expire [29]	This message is displayed by the calibration reminder functionality.	Calibrate the probe and update the calibration date information.
Calibration expired [30]	This warning is displayed by the calibration reminder functionality if the calibration interval has been exceeded.	

More information

- [Technical support \(page 135\)](#)

13.3 Analog output error state

The probe sets the analog output channel into a defined error level instead of the measured result in two situations:

- Probe detects a measurement malfunction. This means an actual measurement problem, such as sensor damage or unsuitable environmental conditions.
- Measured value(s) are significantly outside the scaled output range.

The default error level depends on the output type:

Output	Default error level
0–20 mA	23 mA
4–20 mA	2 mA
0–5 V	0 V
0–10 V	0 V

The probe resumes normal operation of the analog output when the cause of the error state is removed.

14. Technical data

14.1 Product specifications

14.1.1 GMP251 specifications

Table 49 GMP251 measurement performance

Property	Description/value
Measurement range	0–20 %CO ₂
Accuracy ¹⁾	
At 5 %CO ₂	±0.1 %CO ₂
At 0–8 %CO ₂	±0.2 %CO ₂
At 8–20 %CO ₂	±0.4 %CO ₂
Calibration uncertainty	
At 5 %CO ₂	±0.07 %CO ₂
At 20 %CO ₂	±0.27 %CO ₂
Long-term stability	
At 0–8 %CO ₂	±0.3 %CO ₂ /year
At 8–12 %CO ₂	±0.5 %CO ₂ /year
at 12–20 %CO ₂	±1.0 %CO ₂ /year
Temperature dependence	
With compensation at 5 %CO ₂ , 0 ... +50 °C (+32 ... +122 °F)	< ±0.05 %CO ₂
With compensation, 0–20 %CO ₂ , –40 ... +60 °C (–40 ... +140 °F)	±0.045 % of reading/°C
without temperature compensation at 5 %CO ₂ (typical)	–0.25 % of reading/°C
Pressure dependence	
With compensation at 5 %CO ₂ 700–1100 hPa	±0.05 %CO ₂
With compensation, 0–20 %CO ₂ 500–1100 hPa	±0.015 % of reading/hPa
Without compensation (typical)	+0.15 % of reading/hPa
Humidity dependence	

Property	Description/value
With compensation, 0–20 %CO ₂ , 0–100 %RH	±0.7 % of reading (at +25 °C (+77 °F))
Without compensation (typical)	+0.05 % of reading / %RH
O₂ dependence	
With compensation, 0–20 %CO ₂ , 0–90 %O ₂	±0.6 % of reading (at +25 °C (+77 °F))
Without compensation (typical)	-0.08 % of reading / %O ₂
Start-up, warm-up, and response time	
Start-up time at +25 °C (+77 °F)	< 10 s
Warm-up time for full spec.	< 4 min
Response time (T90):	
With standard filter	< 1 min
Flow-through option with > 0.1 l/min	< 1 min
With spray shield	< 2 min
Flow rate dependence (for flow-through option)	
Flow rate dependence:	
< 1 l/min flow	No effect
1–10 l/min flow	< 0.6 % of reading/ l/min

1) At 25 °C (77 °F) and 1013 hPa (incl. repeatability and non-linearity).

Table 50 GMP251 operating environment

Property	Description/value
Operating temperature of CO ₂ measurement	-40 ... +60 °C (-40 ... +140 °F) ¹⁾
Storage temperature	-40 ... +70 °C (-40 ... +158 °F)
Humidity	0–100 %RH, non-condensing
Condensation prevention	Sensor head heating, when power on
IP rating, probe body	IP65: Dust-tight. Protected from water jets from any direction.
Chemical tolerance (temporary exposure during cleaning)	<ul style="list-style-type: none"> • H₂O₂ (2000 ppm, non-condensing) • Alcohol-based cleaning agents (for example ethanol and IPA) • Acetone • Acetic acid
Pressure	

Property	Description/value
Compensated	500–1100 hPa
Operating	< 1.5 bar
Gas flow (for flow-through option)	
Operating range	< 10 l/min
Recommended range	0.1–0.8 l/min

- 1) Occasional short-term exposure to up to +90 °C (+194 °F) allowed, provided that the probe is fully installed inside the measured condition and power is switched on. Accuracy specification not applicable if used in temperatures above +60 °C (+140 °F).

Table 51 GMP251 inputs and outputs

Property	Description/value
Analog outputs	<ul style="list-style-type: none"> 0–5/10 V (scalable), min. load 10 kΩ 0/4–20 mA (scalable), max. load 500 Ω
Digital output	Over RS-485: <ul style="list-style-type: none"> Modbus Vaisala Industrial Protocol
Operating voltage	
With digital output in use	12–30 V DC
With voltage output in use	12–30 V DC
With current output in use	20–30 V DC
Power consumption	
Typical (continuous operation)	0.4 W
Maximum	0.5 W

Table 52 GMP251 mechanical specifications

Property	Description/value
Weight, probe	45 g (1.59 oz)
Connector type	M12 5-pin male
Materials	
Probe housing	PBT polymer
Filter	PTFE membrane, PBT polymer grid
Connector	Nickel plated brass
Dimensions	

Property	Description/value
Probe diameter	25 mm (0.98 in)
Probe length	96 mm (3.78 in)

Table 53 GMP251 compliance

Property	Description/Value
EU directives and regulations	EMC Directive (2014/30/EU) RoHS Directive (2011/65/EU) as amended by 2015/863
Electromagnetic compatibility (EMC)	EN 61326-1, basic electromagnetic environment
EMC emissions	CISPR 32 / EN 55032, Class B
Compliance marks	CE, RCM

14.1.2 GMP80P specifications

Table 54 GMP80P measurement performance

Property	Description/Value
Measurement range	0–20 %CO ₂
Accuracy ¹⁾	
At 5 %CO ₂	±0.1 %CO ₂
At 0–8 %CO ₂	±0.2 %CO ₂
At 8–20 %CO ₂	±0.4 %CO ₂
Calibration uncertainty	
At 5 %CO ₂	±0.07 %CO ₂
At 20 %CO ₂	±0.27 %CO ₂
Long-term stability	
At 0–8 %CO ₂	±0.3 %CO ₂ /year
At 8–12 %CO ₂	±0.5 %CO ₂ /year
At 12–20 %CO ₂	±1.0 %CO ₂ /year
Start-up, warm-up, and response time	
Start-up time at +25 °C (+77 °F)	< 10 s
Warm-up time to full specification	< 4 min

Property	Description/Value
Response time (T90)	< 1 min
Other	
Sampling frequency	1 s

1) At 25 °C (77 °F) and 1013 hPa (incl. repeatability and non-linearity).

Table 55 GMP80P operating environment

Property	Description/Value
Operating environment	Indoor use
Operating temperature	+15 ... +40 °C (+59 ... +104 °F)
Storage temperature	-20 ... +60 °C (-4 ... +140 °F)
Operating humidity	0-95 %RH
Storage humidity	20-85 %RH
IP rating	IP4X
Flow rate of sampled/aspirated gas	0.3 l/min

Table 56 GMP80P powering

Property	Description/Value
Operating voltage	12-30 V DC
Power consumption	Typical: 0.7 W Maximum: 0.9 W
Current consumption	60 mA typical at minimum voltage

Table 57 GMP80P digital communication

Property	Description/Value
Digital output	RS-485

Table 58 GMP80P mechanical specifications

Property	Description/Value
Cable connector type	M12 5-pin A-coded male
Barbed fitting for gas outlet (optional)	For tube with ID 4 mm (approx. 5/32 in)
Dimensions (H × W × D)	

Property	Description/Value
Probe body (incl. cable connectors), without sampling pipe:	114 × 102 × 40 mm (4.49 × 4.02 × 1.57 in)
Depth with sampling pipe:	145 mm (5.71 in)
Weight	410 g (14.46 oz)
Materials	
GMP80P (excluding GMP251 probe):	
Housing	Aluminum and polyamide
Sampling pipe	Stainless steel
Barbed fitting for gas outlet	Nickel-plated brass
Connectors (on the bottom)	Nickel-plated brass
GMP251 probe:	
Housing	PBT polymer
Filter	PTFE membrane, PBT polymer grid
Connector	Nickel-plated brass

Table 59 GMP80P compliance

Property	Description/Value
EU directives and regulations	EMC Directive (2014/30/EU) RoHS Directive (2011/65/EU) as amended by 2015/863 REACH Regulation (EC 1907/2006)
Electromagnetic compatibility (EMC)	EN 61326-1, basic electromagnetic environment FCC part 15 B, Class B ICES-3 / NMB-3 (Class B)
Compliance marks	CE

Table 60 GMP80P manufacturer contact information

Property	Description/Value
Made in Finland by	Vaisala Oyj, Vanha Nurmijärventie 21, FI-01670, Vantaa
UK address	Vaisala Ltd, 6230 Bishops Court, Solihull Parkway, Birmingham Business Park, Birmingham B37 7YB

14.2 Spare parts and accessories



Information on spare parts, accessories, and calibration products is available online at vaisala.com and store.vaisala.com.

14.2.1 GMP251 spare parts and accessories

Table 61 GMP251 spare parts and accessories

Item	Item code
Standard membrane filter	ASM211650SP
Porous sintered PTFE filter	DRW243649SP
Probe connection cable with open wires (1.5 m), shielded	223263SP
Probe connection cable with open wires (1.5 m), shielded	254294SP
Probe connection cable with open wires (3 m), shielded	26719SP
Probe connection cable with open wires (10 m), shielded	216546SP
Probe connection cable with open wires and 90° plug (0.6 m), shielded	244669SP
Probe connection cable with open wires and 90° plug (1.5 m), shielded	255102
MI70 connection cable, M12 5-pin	CBL210472
Connection cable for Indigo80, M12-5F - M12-5M, 1.5 m	272075SP
Flat cable for GMP250 probes, M12 5-pin	CBL210493SP
Indigo USB adapter ¹⁾	USB2
Probe mounting clips (2 pcs)	243257SP
Probe mounting flange	243261SP
Probe holder assembly	ASM213582
Flow-through adapter with gas ports	ASM211697SP
Calibration adapter	DRW244827SP
Spray shield	ASM212017SP

1) Vaisala Insight software for Windows available at www.vaisala.com/insight

14.2.2 GMP80P spare parts and accessories

Table 62 GMP80P spare parts and accessories

Description	Item code / Value
GMP251 CO ₂ probe	Configuration code A1COC0N1 ¹⁾
Cable for probes (M12-M12), 1.5 m (4 ft 11 in)	272075SP
Sampling pipe (stainless steel), 100 mm (3.94 in)	DRW259158SP
Sampling tube (Nafion™ membrane and PUR), 2.1 m (6 ft 11 in) ²⁾	GMTUBESP
Sampling tube (Nafion membrane), 0.6 m (1 ft 12 in)	212807SP
Sampling tube (PUR), 1.5 m (4 ft 11 in)	279471SP
Adapter for sampling tube	279473SP
Barbed fitting for gas outlet	279524SP
Magnetic hanger	ASM214318SP
Indigo USB adapter	USB2

1) See the GMP251 Order Form at docs.vaisala.com for more information.

2) Consists of a 0.6-m (1 ft 12 in) Nafion tube and a 1.5-m (4 ft 11 in) polyurethane tube.

14.3 Dimensions

14.3.1 GMP251 probe dimensions

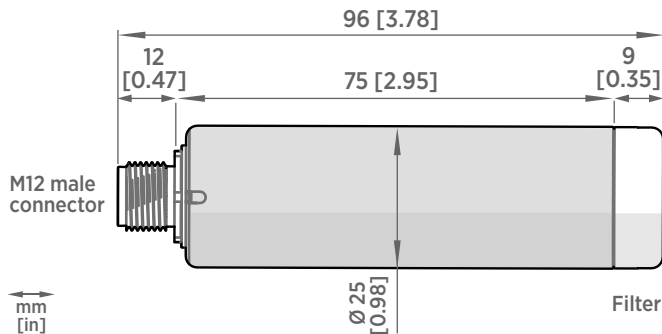


Figure 41 GMP251 dimensions

14.3.2 GMP80P dimensions


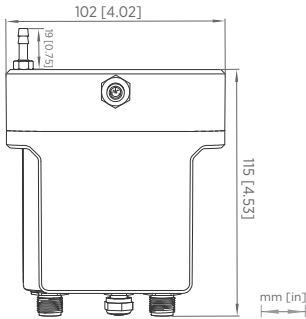


Figure 42 GMP80P dimensions, front and side

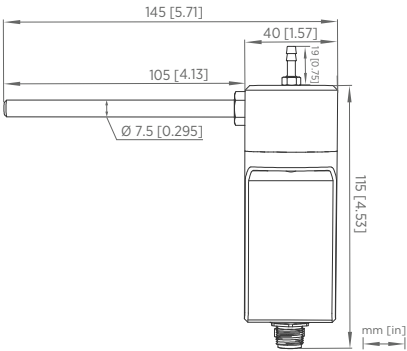

14.3.3 243261SP mounting flange dimensions

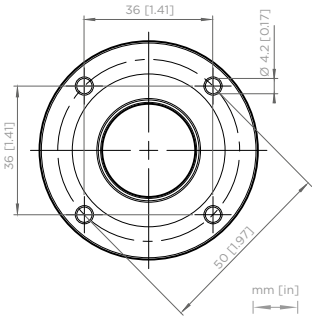


Figure 43 243261SP mounting flange dimensions

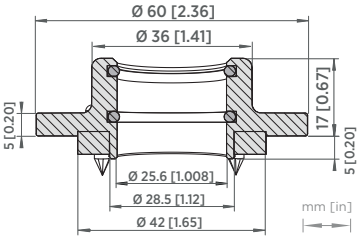


Figure 44 243261SP mounting flange dimensions, cross section

14.3.4 DRW244827SP calibration adapter dimensions

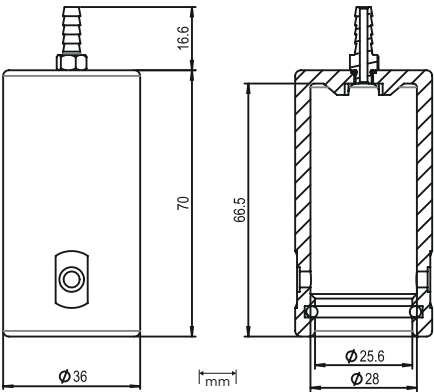


Figure 45 DRW244827SP calibration adapter dimensions

14.3.5 ASM212017SP spray shield dimensions

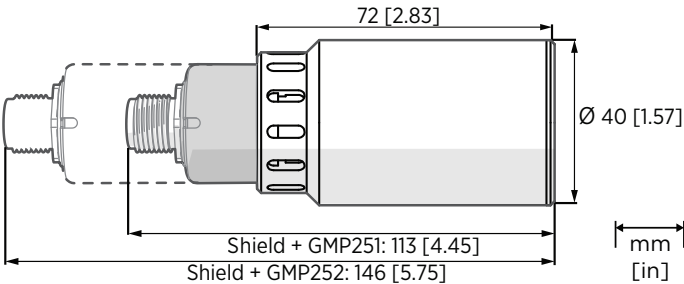


Figure 46 ASM212017SP spray shield dimensions

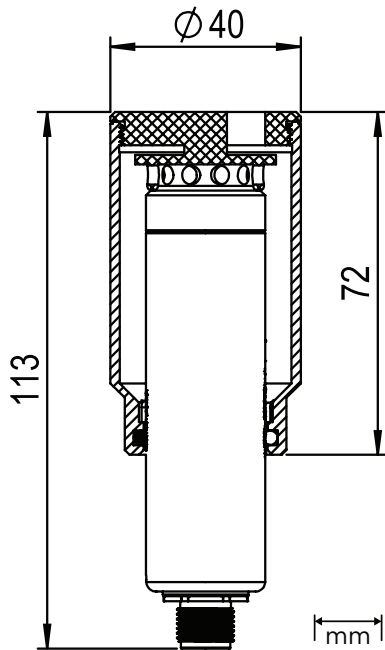


Figure 47 GMP251 spray shield cross section

Appendix A. Modbus reference

A.1 Default communication settings

Table 63 Default Modbus serial communication settings

Property	Description/Value
Serial bit rate	19200
Parity	N
Number of data bits	8
Number of stop bits	2
Modbus device address	240

A.2 Function codes

Table 64 Modbus function codes

Function code (decimal)	Function code (hexadecimal)	Name
03	0x03	Read Holding Registers
16	0x10	Write Multiple Registers
43/14	0x2B / 0x0E	Read Device Identification

A.3 Data encoding

In the data registers, the numeric values are available in one or two formats with separate register addresses: 32-bit IEEE floating point format and/or 16-bit signed integer format.

A.3.1 32-bit floating point or 32-bit integer format

Registers using **32-bit float** data format are encoded using the **binary32** encoding defined in IEEE 754. The format is also known as "single-precision floating point format".

The least significant 16 bits of a floating point number are placed at the Modbus register listed in the table, while the most significant 16 bits are placed in the register with number/address + 1, as specified in Open Modbus TCP Specification, Release 1.0. This is also known as "little-endian" or "Modicon" word order.

Despite the specification, some Modbus masters may expect a "big-endian" word order (most significant word first). In such case, you must select "word-swapped" floating point format in your Modbus master for the Modbus registers of the device.

A complete 32-bit floating point or 32-bit integer value should be read and written in a single Modbus transaction.



CAUTION! Reading the measurement data registers with incorrect floating point format setting may occasionally result in correct-looking, but nevertheless incorrect values.



It is highly recommended to verify that you have configured the floating point format correctly on your Modbus host system by reading a floating point value from a test value register.

A.3.2 16-bit integer format

Some 16-bit integer values in the data registers are scaled to include the necessary decimals. The scaling factors for those values are shown in the register tables.

Table 65 Interpretation of 16-bit signed integer values

Value (decimal)	Value (hexadecimal)	Description
0-32766	0x0000-0x7FFE	Value in range 0-32766
32767	0x7FFF	Value is 32767 or larger
32768	0x8000	Value is not available
32769	0x8001	Value is -32767 or smaller
32770-65535	0x8002-0xFFFF	Value in range -32766 ... -1 (2's complement)

A.4 Modbus registers

Registers are numbered in decimal, starting from 1. Register addresses in actual Modbus messages (Modbus Protocol Data Unit (PDU)) are in hexadecimal and start from zero. Register number 1 corresponds to address 0x0 in the actual Modbus message.



CAUTION! Reading the wrong register(s) may result in correct-looking values. Check the reference documentation of your Modbus host (PLC) to verify which notation it uses for Modbus register addresses.

A.4.1 Measurement data registers

Table 66 Modbus measurement data registers (read-only)

Register number	Address	Description	Data format	Unit
1	0x0000	Measured CO ₂ value	32-bit float	ppm
3	0x0002	Compensation T	32-bit float	°C
5	0x0004	Measured T	32-bit float	°C
257	0x0100	Measured CO ₂ value	16-bit signed integer	ppm (up to 32 000 ppm)
258	0x0101	Measured CO ₂ value	16-bit signed integer	ppm ¹⁾ (scaled, up to approx. 320 000 ppm)

- 1) The ppm output of the second Measured CO₂ value register (number 258) is scaled and must be multiplied by 10.

A.4.2 Configuration registers



CAUTION! Default power-up values (registers 513–519) are written into non-volatile EEPROM memory. The EEPROM memory implementation has a limit of 30000 writes, and is intended to be used only when saving long-term or permanent configurations. Use the volatile memory (registers 521–527, values cleared on power-up) for non-permanent configurations (for example, if the probe is used in a system that regularly updates the compensation values).

Table 67 Modbus configuration data registers (writable)

Register number	Address	Description	Data format	Unit / Valid range
513	0x0200	Power-up value for pressure compensation	32-bit float	hPa 500–1100 hPa Operating <1.5 bar (Init/default: 1013.25)
515	0x0202	Power-up value for temperature compensation	32-bit float	°C –40 ... +80 °C (Init/default: 25)
517	0x0204	Power-up value for humidity compensation	32-bit float	%RH 0–100 %RH (Init/default: 0)

Register number	Address	Description	Data format	Unit / Valid range
519	0x0206	Power-up value for oxygen compensation	32-bit float	%O ₂ 0–100 %O ₂ (Init/default: 0)
521	0x0208	Volatile (value cleared at probe reset) pressure compensation	32-bit float	500–1100 hPa (Init copied from power-up value)
523	0x020A	Volatile (value cleared at probe reset) temperature compensation	32-bit float	–40 ... +80 °C (Init copied from power-up value)
525	0x020C	Volatile (value cleared at probe reset) humidity compensation	32-bit float	0–100 %RH (Init copied from power-up value)
527	0x020E	Volatile (value cleared at probe reset) oxygen compensation	32-bit float	0–100 %O ₂ (Init copied from power-up value)
769	0x0300	Modbus address	16-bit integer	1–247 (default: 240)
770	0x0301	Serial speed	enum	4800–115200 0 = 4800 1 = 9600 2 = 19200 3 = 38400 4 = 57600 5 = 115200 (default: 2 (19200))
771	0x0302	Serial parity	enum	0 = None 1 = Even 2 = Odd (default: 0 (None))
772	0x0303	Serial stop bits	16-bit integer	1–2 (default: 2)
773	0x0304	Pressure compensation mode	enum	0 = Off 1 = On (default: 1 (On))

Register number	Address	Description	Data format	Unit / Valid range
774	0x0305	Temperature compensation mode	enum	0 = Off 1 = Given 2 = Measured (default: 2 (Measured))
775	0x0306	Humidity compensation mode	enum	0 = Off 1 = On (default: 0 (Off))
776	0x0307	Oxygen compensation mode	enum	0 = Off 1 = On (default: 0 (Off))
777	0x0308	CO ₂ filtering factor	16-bit integer	0–100 (default: 100 (no filtering)). For information on setting the filtering factor, see Filtering factor (page 133).



To apply an accurate relative humidity compensation (775), you must also enable temperature compensation (774) and pressure compensation (773).

A.4.3 Status registers

Table 68 Modbus status registers (read-only)

Register number	Address	Description	Data format	Notes
0x2049	0800	Device status	16-bit	0 = Status OK. 1 = Critical error. 2 = Error. 4 = Warning.
0x2050	0801	CO ₂ status	16-bit	0 = Status OK. 2 = CO ₂ reading not reliable. Appears during transmitter start-up. 256 = Measurement not ready. Appears during transmitter start-up.

If several errors are active at the same time, the Modbus error code register contains the sum of currently active error codes.

More information

- [Error messages \(page 111\)](#)

A.5 Device identification objects

Table 69 Device identification objects

Object ID	Object ID (hexadecimal)	Object name	Example contents
0	0x00	VendorName	"Vaisala"
1	0x01	ProductCode	GMP251
2	0x02	MajorMinorVersion	"1.2.3" Software version of the device.
3	0x03	VendorUrl	" http://www.vaisala.com/ "
4	0x04	ProductName	"GMP251 Carbon Dioxide Probe"
128	0x80	SerialNumber ¹⁾	"K0710040"
129	0x81	CalibrationDate ¹⁾	"2020-01-31" Calibration date in YYYY-MM-DD format. Empty string if not set/valid.
130	0x82	CalibrationText ¹⁾	"Vaisala/HEL" Calibration information text. Empty string if not set/valid.

1) *Vaisala-specific device information.*

A.6 Exception responses

Table 70 Modbus exception responses

Code	Name	Reason
01	ILLEGAL FUNCTION	Unsupported function code
02	ILLEGAL DATA ADDRESS	Register address or Object ID out of valid ranges
03	ILLEGAL DATA VALUE	Otherwise invalid request

Accessing unavailable (temporarily missing) measurement data does not generate a Modbus exception. “Unavailable” value (a quiet NaN for floating point data or 0x8000 for integer data) is returned instead. An exception is generated only for any access outside the applicable register ranges.

A.7 Modbus communication examples

Reading CO₂ value




Device address used in the following examples is 240 (0xF0). The values returned by the device differ depending on the ambient conditions and/or device settings. Your device might not return exactly same values.

Request		Response	
Bytes on the line (hexadecimal)	Description	Bytes on the line (hexadecimal)	Description
(silence for 3.5 bytes)	Start of Modbus RTU frame	(silence for 3.5 bytes)	Start of Modbus RTU frame
0xF0	GMP25x address	0xF0	GMP25x address
0x03	Function (Read Holding Registers)	0x03	Function (Read Holding Registers)
0x00	Register address	0x04	Number of data bytes
0x00		0xD4	Value of first register (least significant word)
0x00		0x7A	
0x02	Number of 16-bit registers to read (2)	0x43	Value of second register (most significant word)
0xD1	Modbus RTU checksum	0xE8	Modbus RTU checksum
0x2A		0x33	
(silence for 3.5 bytes)	End of Modbus RTU frame	0xAB	Modbus RTU checksum
		(silence for 3.5 bytes)	End of Modbus RTU frame

Communication description	
Register address	1 (1-based Modbus documentation format) = 0x0000 (0-based format used in actual communication).

Communication description	
Data format	Two 16-bit Modbus registers interpreted as IEEE 754 binary32 floating point value, least significant word first.
Returned value	43E8D47Ah, which is binary32 representation of 465.65997 (ppm).

Writing volatile compensation pressure value

Request		Response		
Bytes on the line (hexadecimal)	Description	Bytes on the line (hexadecimal)	Description	
(silence for 3.5 bytes)	Start of Modbus RTU frame	(silence for 3.5 bytes)	Start of Modbus RTU frame	
0xF0	GMP25x address	0xF0	GMP25x address	
0x10	Function (Write Multiple Registers)	0x10	Function (Write Multiple Registers)	
0x02	Register address	0x02	Register address	
0x08		0x08		
0x00	Number of registers to write (2)	0x00	Number of 16-bit registers written (2)	
0x02		0x02		
0x04	Number of data bytes	0xD4	Modbus RTU checksum	
0x50	Value for the first register (least significant word)	0x93		
0x00		(silence for 3.5 bytes)	End of Modbus RTU frame	
0x44	Value for the second register (most significant word)	<div><div></div><div><p>The response to a write function informs that the function was correctly received by the device. It does not guarantee that the written value was accepted by the device (for example, in case of out-of-range values).</p><p>To verify that the value was really accepted by the device, read the register value after writing.</p></div></div>		
0x7D				
0x0E	Modbus RTU checksum			
0xB7				
(silence for 3.5 bytes)	End of Modbus RTU frame			

Communication description	
Register address	521 (1-based Modbus documentation format) = 0x0208 (0-based format used in actual communication).
Data format	Two 16-bit Modbus registers interpreted as IEEE 754 binary32 floating point value, least significant word first.
Value to write	1013.25 (hPa), in binary32 format 0x447D5000.

A.8 Filtering factor

Modbus register 777 sets the CO₂ filtering factor.

The filtering factor affects the speed at which the latest CO₂ measurement is integrated into the output of the probe. A new measurement is produced approximately every two seconds.

By default, the filtering factor is set to 1.0, which means the latest measurement is shown directly in the output, without any filtering. If the measuring environment produces occasional exceptionally high or low readings that need to be averaged out in the output, filtering can be applied.

To apply filtering, you need to set a filtering factor that determines how much the previous measurements affect the calculation of measurement output. For example, when using filtering factor of 0.1, the new output is a combination of previous measurements (90 %) and the latest measurement (10 %).

Examples of the effect of filtering on output:

- Filtering factor 1.0 = No filtering, the latest measurement is output directly without integrating previous measurements.
- Filtering factor 0.5 = The reading output shows ~75 % of the measurement change after two 2-second measurement cycles and ~93 % after four cycles.
- Filtering factor 0.1 = The reading output shows ~90 % of the measurement change after 22 measurement cycles.

The configuration range of the filtering factor is 0–100 in the 16-bit register: for example, to set the factor to 0.5, set the value of the register to 50.

The following formula is used when calculating the output:

$$o_{\text{new}} = o_{\text{old}} + (m_{\text{new}} - o_{\text{old}}) \times f$$

o_{new}	New output
o_{old}	Previous output
m_{new}	New measurement
f	Filtering factor

Appendix B. CO₂ dilution coefficients reference

B.1 CO₂ dilution coefficients

Table 71 (page 134) lists the dilution coefficients for CO₂ concentration when drying a gas sample using Nafion tubing. The table lists the coefficients for dew point temperature scale -40 ... +60 °C. A similar table for a narrower scale (0 ... +60 °C) is presented in [CO₂ dilution coefficients for pump sampling method \(page 54\)](#).

Table 71 CO₂ dilution coefficients for measurements with Nafion tube, T_d of sampled gas -40 ... +60 °C

	Dew point temperature of sampled gas (°C)										
T _d (°C) ¹⁾	-40	-30	-20	-10	0	+10	+20	+30	+40	+50	+60
-60	0.9999	0.9996	0.999	0.997	0.994	0.998	0.997	0.958	0.927	0.878	0.803
-50	0.9999	0.9997	0.999	0.997	0.994	0.998	0.997	0.958	0.927	0.878	0.803
-40	1.000	0.9998	0.999	0.998	0.994	0.998	0.997	0.958	0.927	0.878	0.803
-30		1.000	0.999	0.998	0.994	0.998	0.997	0.958	0.928	0.879	0.804
-20			1.000	0.998	0.995	0.999	0.978	0.959	0.928	0.879	0.804
-10				1.000	0.997	0.990	0.979	0.961	0.930	0.881	0.806
0					1.000	0.994	0.983	0.964	0.933	0.884	0.809
+10						1.000	0.989	0.970	0.939	0.890	0.815
+20							1.000	0.981	0.950	0.901	0.826
+30								1.000	0.969	0.920	0.845
+40									1.000	0.951	0.876
+50										1.000	0.925
+60											1.000

1) Dew point temperature of environment where GMP80P is located.

Warranty

For standard warranty terms and conditions, see vaisala.com/warranty.

Please observe that any such warranty may not be valid in case of damage due to normal wear and tear, exceptional operating conditions, negligent handling or installation, or unauthorized modifications. Please see the applicable supply contract or Conditions of Sale for details of the warranty for each product.

Technical support



Contact Vaisala technical support at helpdesk@vaisala.com. Provide at least the following supporting information as applicable:

- Product name, model, and serial number
- Software/Firmware version
- Name and location of the installation site
- Name and contact information of a technical person who can provide further information on the problem

For more information, see vaisala.com/support.

Recycling



Recycle all applicable material according to local regulations.

VAISALA

